



(11) **EP 1 425 294 B1**

(12) **EUROPEAN PATENT SPECIFICATION**

(45) Date of publication and mention
of the grant of the patent:
16.07.2008 Bulletin 2008/29

(51) Int Cl.:
G01N 33/543 (2006.01) G01N 33/569 (2006.01)

(21) Application number: **02761478.3**

(86) International application number:
PCT/US2002/026861

(22) Date of filing: **23.08.2002**

(87) International publication number:
WO 2003/019141 (06.03.2003 Gazette 2003/10)

(54) **ANALYSIS OF CIRCULATING TUMOR CELLS, FRAGMENTS, AND DEBRIS**

ANALYSE ZIRKULIERENDER TUMORZELLEN, FRAGMENTE UND TRÜMMER

ANALYSE DE CELLULES TUMORALES EN CIRCULATION, DE FRAGMENTS ET DE DEBRIS
ASSOCIES

(84) Designated Contracting States:
**AT BE BG CH CY CZ DE DK EE ES FI FR GB GR
IE IT LI LU MC NL PT SE SK TR**

(30) Priority: **23.08.2001 US 314151 P**
03.04.2002 US 369628 P

(43) Date of publication of application:
09.06.2004 Bulletin 2004/24

(73) Proprietor: **IMMUNIVEST CORPORATION**
Wilmington, DE 19899 (US)

(72) Inventors:
• **RAO, Galla, Chandra**
Princeton Junction, NJ 08550 (US)
• **LARSON, Christopher**
Media, PA 19063 (US)
• **REPOLLET, Madeline**
Ft. Washington, PA 19034 (US)
• **RUTNER, Herman**
Hatboro, PA 19040 (US)
• **TERSTAPPEN, Leon**
Huntingdon Valley, PA 19006 (US)
• **O'HARA, Shawn, Mark**
Ambler, PA 19002 (US)

• **GROSS, Steven**
Ambler, PA 19002 (US)

(74) Representative: **Hanson, William Bennett et al**
Bromhead Johnson
19 Buckingham Street
London WC2N 6EF (GB)

(56) References cited:
WO-A-00/47998 WO-A-01/19956
WO-A-99/41613 US-A- 5 385 707
US-A1- 2001 024 802 US-B1- 6 190 870
US-B1- 6 365 362

• **BRANDT ET AL.: 'Isolation of prostate-derived
single cells and cell clusters from human
peripheral blood' CANCER RESEARCH vol. 56,
no. 20, 15 October 1996, pages 4556 - 4561,
XP001155823**
• **TANAKA ET AL.: 'Immunohistochemical study of
localization of extracellular matrix after holmium
YAG laser irradiation rat cornea' JAPANESE
JOURNAL OF OPHTHALMOLOGY vol. 44, no. 5,
September 2000 - October 2000, pages 482 - 484,
XP002972404**

Note: Within nine months of the publication of the mention of the grant of the European patent in the European Patent Bulletin, any person may give notice to the European Patent Office of opposition to that patent, in accordance with the Implementing Regulations. Notice of opposition shall not be deemed to have been filed until the opposition fee has been paid. (Art. 99(1) European Patent Convention).

EP 1 425 294 B1

Description

BACKGROUND OF THE INVENTION

5 **[0001]** Many clinicians believe that cancer is an organ-confined disease in its early stages. However, it appears that this notion is incorrect, and cancer is often a systemic disease by the time it is first detected using methods currently available. There is evidence that primary cancers begin shedding neoplastic cells into the circulation at an early disease stage prior to the appearance of clinical manifestations. Upon vascularization of a tumor, tumor cells shed into the circulation may attach and colonize at distant sites to form metastases. These circulating tumor cells (CTC) contain
10 markers not normally found in healthy individuals' cells, thus forming the basis for diagnosis and treatment of specific carcinomas. Hence, the presence of tumor cells in the circulation can be used to screen for cancer in place of, or in conjunction with, other tests, such as mammography for breast cancer, or measurements of Prostate Specific Antigen (PSA) for prostate cancer. By employing appropriate monoclonal antibodies directed to associated markers on or in target cells, or by using other assays for cell protein expression, or by the analysis of cellular mRNA, the organ origin of such cells may readily be determined, e.g., breast, prostate, colon, lung, ovarian or other non-hematopoietic cancers.

15 **[0002]** Thus, in cases where cancer cells can be detected, while there are essentially no clinical signs of a tumor, it will be possible to identify their presence as well as the organ of origin. If CTC are detected after surgery, one could stat adjuvant therapy if the CTC are a sign of relapse. Predicting the patient's need for such treatment, or the efficacy thereof, given the costs of such therapies, is a significant and beneficial piece of clinical information. It is also clear that if the number of CTC changes, it may predict progression (CTC increase) or response (CTC decrease).
20 **[0003]** Malignant tumors are characterized by their ability to invade adjacent tissue. In general, tumors with a diameter of 1mm are vascularized and animal studies show that as much as 4% of the cells present in the tumor can be shed into the circulation in a 24 hour period (Butler, TP & Gullino PM, 1975 Cancer Research 35:512-516). The shedding capacity of a tumor is most likely dependent on the aggressiveness of the tumor. Although tumor cells are shed into the circulation on a continuous basis, it is believed that none or only a small fraction will give rise to distant metastasis (Butler & Gullino, *supra*). Increase in tumor mass might be expected to be proportional to an increase in the frequency of CTC, as single cells, or clusters of cells with increased adhesiveness (and possibly greater invasive potential). If this were found to be the case, methods available with a high level of sensitivity would facilitate assessment of tumor load in patients with distant metastasis as well as those with localized disease. Detection of tumor cells in peripheral blood of patients with localized disease has the potential not only to detect a tumor at an earlier stage but also to provide indications as to the potential invasiveness of the tumor.

25 **[0004]** However, whole blood is a complex body fluid containing diverse populations of cellular and soluble components capable of undergoing numerous biochemical and enzymatic reactions *in vivo* and *in vitro*, particularly on prolonged storage for more than 24hrs. Some of these reactions are related to immunological destruction of CTC that are recognized as foreign species. The patient's immune response weakens or destroys tumor cells by the normal defense mechanisms including phagocytosis and neutrophil activation.
30 Chemotherapy similarly is intended to reduce both cell function and proliferation by inducing cell death by necrosis. Besides these external destructive factors, tumor cells damaged in a hostile environment may undergo programmed death or apoptosis. Normal and abnormal cells (including CTC) undergoing apoptosis or necrosis, have altered membrane permeabilities that allow escape of DNA, RNA, and other intracellular components leading to formation of damaged cells, fragmented cells, cellular debris, and eventual complete disintegration. Such tumor cell debris may still bear epitopes or determinants characteristic of intact cells and can lead to spurious increases in the number of detected circulating cancer cells. Whole blood specimens from healthy individuals also have been observed to undergo destruction of labile blood cell components, herein categorized as decreased blood quality, on prolonged storage for periods of greater than 24 hours. For example, erythrocytes may rupture and release hemoglobin and produce cell ghosts. Leukocytes, particularly granulocytes, are known to be labile and diminish on storage. Such changes increase the amount of cellular debris that can interfere with the isolation and detection of rare target cells such as CTC. The combined effects of these destructive processes, collectively defined as post-draw or *in vitro* damage, can substantially increase cellular debris, which is readily detectable, for instance, in flow cytometric and microscopic analyses, such as CellSpotter® or CellTracks™.

35 **[0005]** Detection of circulating tumor cells by microscopic imaging is similarly adversely affected by spurious decreases in classifiable tumor cells and a corresponding increase in interfering stainable debris. Hence, maintaining the integrity or the quality of the blood specimen is of utmost importance, since there may be a delay of as much as 24 hours between blood draw and specimen processing. Such delays are to be expected, since the techniques and equipment used in processing blood for this assay may not be readily available in every laboratory. The time necessary for a sample to arrive at a laboratory for sample processing may vary considerably. It is therefore important to establish the time window within which a sample can be processed. In routine hematology analyses, blood samples can be analyzed within 24 hours. However, as the analysis of rare blood cells is more critical, the time window in which a blood sample can be
40
45
50
55

EP 1 425 294 B1

analyzed is shorter. An example is immunophenotyping of blood cells, which, in general, must be performed within 24 hours. In a cancer cell assay, larger volumes of blood have to be processed, and *in vitro* degradation of the blood sample can become more problematic as materials released by disintegrating cells, both from CTC and from hematopoietic cells, can increase the background and therefore decrease the ability to detect tumor cells.

5 **[0006]** The origin and nature of observed small debris and large clump-like aggregates are not fully understood, but are believed to involve cellular components or elements originating from target cells, non-target cells, and possibly plasma components. Since CTC can be considered immunologically foreign species, normal cellular immune responses of the host will occur *in vivo* even before blood draw. Also large numbers of CTC can be continuously shed from a tumor site, and a steady-state level is maintained in which destruction of CTC equals the shedding rate which in turn depends on the size of the tumor burden (see JG Moreno *et al.* "Changes in Circulating Carcinoma Cells in Patients with Metastatic Prostate Cancer Correlates with Disease State." *Urology* **58**, 2001).

10 **[0007]** Various methods are known in this particular art field for recovering tumor cells from blood. For example, US Patent #6,190,870 to AmCell and Miltenyi teaches immunomagnetic isolation followed by flow cytometric enumeration. However, before immunomagnetic separation, the blood samples are pre-processed using density gradients. Further-
15 more, there is no discussion of isolating or counting anything other than intact cells. There is also no visual morphological analysis of the samples.

[0008] In US Patent #6,197,523, Rimm *et al.* describe enumerating cancer cells in 100 μ l blood samples. The methods use capillary microscopy to confirm the identity of cells that are found. The methods are specific for intact cells, which must be present in high numbers due to small sample volume. However, there is no discussion of isolating or enumerating
20 anything other else, such as fragments or debris.

[0009] In US Patent #6,365,362 to Immunivest, methods are described for immunomagnetically enriching and analyzing samples for tumor cells in blood. The methods are specifically directed towards analyzing intact cells, where the number of cells correlates with the disease state. The isolated cells are labeled for the presence of nucleic acid and an additional marker, which allows the exclusion of non-target sample components during analysis.

25 **[0010]** In WO02/20825, Chen describes the use of an adhesion matrix for enumerating tumor cells. Briefly, the matrix is coated with specific adhesion molecules that will bind viable cancer cells with presumed metastatic potential. The matrix can then be analyzed for the presence and type of captured cells. Also described are methods for using the matrix in screening treatments. While steps are taken to discriminate between intact cells and apoptotic or necrotic cells, the apoptotic or necrotic cells are specifically excluded from analysis.

30 **[0011]** In WO00/47998 from Cell Works, two pathways are described for CTC, terminal and proliferative. Both pathways begin with an "indeterminate" cell that progresses, as determined by morphological differences, down either the terminal or proliferative pathway. A cell in the terminal pathway eventually is destroyed, and a cell in the proliferative pathway will form a new metastatic colony as a metastatic tumor. This proliferative pathway (cells and clusters) is the focus as a potential diagnostic, but minimizes the significance of cell fragments or debris. These two pathways were designed to
35 explain morphological differences seen in patient samples.

[0012] Generally, the more resistant and proliferative cells survive to establish secondary or metastatic sites. In the peripheral circulation, CTC are further attacked *in vivo* (and also *in vitro*) by activated neutrophils and macrophages resulting progressively in membrane perforation, leakage of electrolytes, smaller molecules, and eventual loss of critical cellular elements including DNA, chromatin, etc, which are essential for cell viability. At a critical point of the cell's demise,
40 cell destruction is further assisted by apoptosis. Apoptosis is characterized by a series of stepwise slow intracellular events, which differs from necrosis or rapid cell death triggered or mediated by an extracellular species, e.g. a cytotoxic anti-tumor drug. All or some of these destructive processes may lead to formation of debris and/or aggregates including stainable DNA, DNA fragments and "DNA ladder" structures from disintegrating CTC as well as from inadvertent destruction of normal hematopoietic cells during drug therapy, since most cytotoxic drugs are administered at near toxic
45 doses.

[0013] As shown in WO00/47998, US #6,190,870, and other publications, CTC can circulate as both live and dead cells, wherein "dead" comprises the full range of damaged and fragmented cells as well as CTC-derived debris. The tumor burden is probably best represented by the total of both intact CTC, including clusters, and damaged CTC, which bear morphological characteristics of cells, but are distinct from clumps and/or aggregates. However, some damaged
50 cells, may have large pores allowing leakage of the liquid and particulate cytosolic contents resulting in a change in the buoyant densities from about 1.06-1.08 to greater than 1.12, or well above the densities of RBC (live and dead cells can be separated at the interface of gradients of $d=1.12$ and 1.16 according to a Pharmacia protocol). Conventional density gradients, as used in # WO00/479 would lose such damaged CTC in the discarded RBC layer having a range in density of about 1.08 to 1.11. CTC debris that is positively stained for cytokeratin may also have densities falling in the RBC or higher ranges, since most intracellular components (with the possible exception of lipophilic membrane fragments that may be located near the plasma-buffy coat interface) have densities in the range of 1.15 to 1.3. Hence, a substantial portion of damaged CTC and CTC debris may be located outside the buffy coat layer, and would not be seen by the density gradient methods, such as those in WO00/47998. Some images of damaged or fragmented CTC are shown,
55

EP 1 425 294 B1

but it is quite possible the damage occurred during cytopspin or subsequent processing, and is thus artifactual. While the densities of most intact tumor cells may fall in the WBC region, it is quite likely that damaged CTC in patient samples have higher densities that may place them in the RBC layer; outside the reach of gradient techniques.

[0014] US Patent Application #2001/0024802 describes methods for binding fragments and debris to beads. That published application described numerous possibilities for the density of fragments and debris of interest. Upon centrifugation, the beads will be located in a layer above RBC, because of the pre-determined specific gravity (density) of the beads coupled to fragments and/or debris. However, this system is dependent on correctly binding fragments and debris to these beads. If any other sample component binds the beads, they may not appear in the desired location, and subsequently will not be subject to analysis.

[0015] Epithelial cells in their tissue of origin obey established growth and development "rules". Those rules include population control. This means that under normal circumstances the number and size of the cells remains constant and changes only when necessary for normal growth and development of the organism. Only the basal cells of the epithelium or immortal cells will divide and they will do so when it is necessary for the epithelium to perform its function, whatever it is depending in the nature and location of the epithelium. Under some abnormal but benign circumstances, cells will proliferate and the basal layer will divide more than usual, causing hyperplasia. Under some other abnormal but benign circumstances, cells may increase in size beyond what is normal for the particular tissue, causing cell gigantism, as in folic acid deficiency.

[0016] Epithelial tissue may increase in size or number of cells also due to pre-malignant or malignant lesions. In these cases, changes similar to those described above are accompanied by nuclear abnormalities ranging from mild in low-grade intraepithelial lesions to severe in malignancies. It is believed that changes in these cells may affect portions of the thickness of the epithelium and as they increase in severity will comprise a thicker portion of such epithelium. These cells do not obey restrictions of contact inhibition and continue growing without tissue controls. When the entire thickness of the epithelium is affected by malignant changes, the condition is recognized as a carcinoma *in situ* (CIS).

[0017] The malignant cells eventually are able to pass through the basement membrane and invade the stroma of the organ as their malignant potential increases. After invading the stroma, these cells are believed to have the potential for reaching the blood vessels. Once they infiltrate the blood vessels, cells find themselves in a completely different environment from the one they originated from.

[0018] The cells may infiltrate the blood vessels as single cells or as clusters of two or more cells. A single cell of epithelial origin circulating through the circulatory system is destined to have one of two outcomes. It may die or it may survive.

Single Cells:

1. The cell may die either through apoptosis due to internal changes or messages in the cell itself. These messages may have been in the cell before intravasation or they may be received while in the blood, or it may die due to the influence of the immune system of the host, which may recognize these cells as "alien" to this environment. The results of cellular death are identifiable in CellSpotter® as enucleated cells, speckled cells or amorphous cells. These cells do not have the potential for cell division or for establishing colonies or metastases.

- Enucleated cells are the result of nuclear disintegration and elimination (karyorrhexis and karyolysis). They are positive for cytokeratin, and negative for nucleic acid.
- The speckled cells are positive for cytokeratin and DAPI and show evidence of cellular degeneration and cytoplasmic disintegration. These cells may represent response to therapy or to the host's immune system as the cytoskeletal proteins retract.
- Another dying tumor cell identifiable using CellSpotter® is the amorphous cell. These cells are probably damaged during the preparation process, a sign that these may be weaker, more delicate cells but may also be the result of apoptosis or immune attack.

2. A viable malignant epithelial cell may have the potential to survive the circulation and form colonies in distant organs. These "survivor cells" appear in CellSpotter® as intact cells with high nuclear material/cytoplasmic material ratio. These cells are probably undifferentiated and can potentially divide in blood and form small clusters that may extravasate in a distant capillary, where the cell may establish a new colony, or it may remain as a single cell until it extravasates, dividing once it establishes itself in the new tissue, starting this way a new colony.

Clusters: The primary tumor may shed clusters that enter the circulation as described by B Brandt et al. ("Isolation of prostate-derived single cells and cell clusters from human peripheral blood." Cancer Research 56 p4556-4561. 1996). These clusters may remain as clusters and invade a distant tissue or they may become dissociated in the circulation, probably due to differences in pressure in blood or to the immune system's intervention. If these cells are dissociated into single cells, they may follow one of the two paths described for single cells above (see 1 and 2). Cluster formations

EP 1 425 294 B1

may have an effect in survival by using the outside cells as a shield that protects the inner cells from the immune system.

[0019] Once a new colony is established in a new organ, some malignant cells will continue replicating to form a new tumor. If they reach new capillaries, the metastasis story may be repeated and secondary metastasis occurs.

Monitoring of treatment in patients with known carcinomas: A decrease in the number of tumor cells and/ or increase in the response index may represent a response to patient therapy.

- Total tumor cells = Dying cells + Survivor cells ($TTC = DC + SC$)
- Response Index = dying cells / total tumor cells ($RI = DC / TTC$).

[0020] The higher the response index, the better the response to therapy. A low response index may indicate that the patient is not responding to the treatment and or that the pt's immune system is not able to handle the tumor load.

[0021] A patient who has 50 total tumor cells that were all survivor cells at pre-treatment visit ($a\ RI = 0/50 = 0$) and has 50 TTC on follow-up (after treatment) visit may have different outcomes depending in the RI. If all the TTC are SC (i.e. $DC = 0$), there was no response to therapy. If there are 50 cells but the response index is $40/50 = 0.8$, then either the immune system or the therapy is having a negative effect on tumor load, therefore, is a positive patient response.

Decisions in follow-up on patients with known pre-malignancies: When a pap smear is diagnosed as having cells with atypia or low-grade intraepithelial lesions, there is always the possibility that these patients have a more severe abnormality, which cells were missed as a sampling error. These patients can be colposcoped and biopsied or they may be asked to return in three months for a repeat pap smear. If the atypical cells were concurrent with a small focal area of malignant cells that did not get sampled, the patient will wait 3 months before she gets any follow-up. This may explain why some pre-malignancies seem to progress quicker than others (misdiagnoses due to sampling error, causing an artifact in statistics). These are usually explained as being a more "aggressive" pre-malignancy. CellSpotter® can be used to help in the decision tree of these patients. All patients with an abnormal pap (5-10% of the pap smears in the USA) can immediately be tested for circulating epithelial cells. Patients with positive tests should be followed-up immediately and aggressively. Patients with negative results may wait the three months for the repeat pap. This would simplify the decision making process for the physician and health professionals and help the patient trust her follow-up procedure. Screening: CellSpotter® image analysis may be used for screening of the general population with the condition that special, tissue specific antibodies would be used on a second test on all abnormal samples. Identification of CTC in a patient may indicate that there is a primary malignancy that has started or is starting the process of metastasis. If these cells are identified as of the tissue of origin with new markers, then organ specific tests, like CT guided fine needle aspirations (FNA) can be used to verify the presence or absence of such malignancies. Patients where a primary cannot be identified may be followed-up with repeat tests after establishing an individual base line.

[0022] In summary, all or some of the above-cited factors can and were found to contribute to debris and/or aggregate formation that have been observed to confound the detection of CTC by direct enrichment procedures from whole blood as disclosed in this invention. The number of intact CTC, damaged or Suspect CTC as well as the degree of damage to the CTC, may further serve as diagnostically important indicators of the tumor burden, the proliferative potential of the tumor cells and/or the effectiveness of therapy. In contrast, the methods and protocols of the prior art combine unavoidable *in vivo* damage to CTC with avoidable *in vitro* storage and processing damage, thus yielding erroneous information on CTC and tumor burdens in cancer patients. Finally, the relatively simple blood test of the present invention described herein, which functions with a high degree of sensitivity and specificity, can be thought of as a "whole body biopsy."

[0023] WO 99/41613 describes methods and reagents for detecting, enumerating and characterizing intact rare carcinoma cells in the blood.

BRIEF DESCRIPTION OF THE INVENTION

[0024] The methods and reagents described in this invention are used to analyze circulating tumor cells, clusters, fragments, and debris- Analysis is performed with a number of platforms, including flow cytometry and the CellSpotter® fluorescent microscopy imaging system. The Examples show the importance of not only analyzing Obvious or Intact CTC, but Suspect CTC or damaged fragments, clusters of CTC, and debris. It is possible to mimic the damage that forms fragments and debris. It is also possible to inhibit further damage of CTC between the blood draw and sample processing through the use of stabilizing agents.

[0025] It has been shown herein that the ability to differentiate between *in vitro* damage, caused by specimen acquisition, transport, storage, processing, or analysis, and *in vivo* damage, caused by apoptosis, necrosis, or the patient's immune system. Indeed, it is desirable to confine, reduce, eliminate, or at least qualify *in vitro* damage to prevent it from interfering in analysis.

[0026] The invention provides a method of testing for the presence of CTCs according to claim 1 and a method of monitoring the presence of CTCs according to claim 2.

EP 1 425 294 B1

BRIEF DESCRIPTION OF THE FIGURES

[0027]

Figure 1 - Models of tumor shedding and metastasis. 1a, shows possible stages of cells, clusters, and fragments. 1b. shows the same model with actual images from samples.

Figure 2-Flow cytometric analysis of immunomagnetically enriched tumor cells from 7.5ml blood.

Figure 3 - CellSpotter® analysis of a 7.5ml blood sample from a metastatic prostate cancer patient that was immunomagnetically enriched for tumor cells. The lines of thumbnails correspond to the different dyes used in the staining process showing tumor candidates stained with cytokeratin PE (green) and DAPI (magenta).

Figure 4 - CellSpotter® classifications of tumor cells isolated from a single whole blood sample of a patient with metastatic prostate cancer stained with cytokeratin PE (green) and DAPI (magenta).

- A - intact cells
- B - damaged tumor cells
- C - tumor cell fragments

Figure 5 - A comparison of the number of Obvious CTC to Suspect CTC in 20 clinical samples.

Figure 6 - CellSpotter® classifications of paclitaxel treated LnCaP cells spiked into whole blood and isolated then stained with cytokeratin PE (green) and DAPI (magenta).

- A - intact cells
- B - dying tumor cells
- C - tumor cell fragments

DETAILED DESCRIPTION OF THE INVENTION

[0028] Herein, various terms that are well understood by those of ordinary skill in the art are used. The intended meaning of these terms does not depart from the accepted meaning.

[0029] The evidence that minimal residual disease in patients with carcinoma has clinical significance is mounting. To effectively monitor minimal residual disease, a qualitative and quantitative assessment is needed. As the frequency of carcinoma cells in blood or bone marrow is low, the laborious manual sample preparation methods involved in the preparation of samples for analysis often leads to erroneous results. To overcome these limitations a semi-automated sample preparation system was developed that minimize variability and provide more consistent results, as described in commonly-owned pending US Application No. 10/081,996, filed 20 February 2002.

[0030] Various methods are available for analyzing or separating the above-mentioned target substances based upon complex formation between the substance of interest and another substance to which the target substance specifically binds. Separation of complexes from unbound material may be accomplished gravitationally, e.g. by settling, or, by centrifugation of finely divided particles or beads coupled to the target substance. Such particles or beads may be made magnetic to facilitate the bound/free separation step. Magnetic particles are well known in the art, as is their use in immune and other bio-specific affinity reactions. Generally, any material that facilitates magnetic or gravitational separation may be employed for this purpose. However, it has become clear that magnetic separation means are the method of choice.

[0031] Magnetic particles can be classified on the basis of size; large (1.5 to about 50 microns), small (0.7-1.5 microns), or colloidal (<200nm), which are also referred to as nanoparticles. Nanoparticles, also known as ferrofluids or ferrofluid-like materials, have many of the properties of classical ferrofluids, and are sometimes referred to herein as colloidal, superparamagnetic particles.

[0032] Small magnetic particles of the type described above are quite useful in analyses involving bio-specific affinity reactions, as they are conveniently coated with biofunctional polymers (e.g., proteins), provide very high surface areas and give reasonable reaction kinetics. Magnetic particles ranging from 0.7-1.5 microns have been described in the patent literature, including, by way of example, US Patent Nos. 3,970,518; 4,018,886; 4,230,685; 4,267,234; 4,452,773; 4,554,088; and 4,659,678. Certain of these particles are disclosed to be useful solid supports for immunological reagents.

[0033] The efficiency with which magnetic separations can be done and the recovery and purity of magnetically labeled cells will depend on many factors. These include:

- number of cells being separated,
- receptor or epitope density of such cells,
- magnetic load per cell,

EP 1 425 294 B1

- non-specific binding (NSB) of the magnetic material,
- carry-over of entrapped non-target cells,
- technique employed,
- nature of the vessel,
- 5 • nature of the vessel surface,
- viscosity of the medium, and
- magnetic separation device employed.

10 **[0034]** If the level of non-specific binding of a system is substantially constant, as is usually the case, then as the target population decreases so will the purity.

[0035] Less obvious is the fact that the smaller the population of a targeted cell, the more difficult it will be to magnetically label and to recover. Furthermore, labeling and recovery will markedly depend on the nature of magnetic particle employed. For example, when cells are incubated with large magnetic particles, such as Dynal beads, cells are labeled through collisions created by mixing of the system, as the beads are too large to diffuse effectively. Thus, if a cell were
15 present in a population at a frequency of 1 cell per ml of blood or even less, as may be the case for tumor cells in very early cancers, then the probability of labeling target cells will be related to the number of magnetic particles added to the system and the length of time of mixing. Since mixing of cells with such particles for substantial periods of time would be deleterious, it becomes necessary to increase particle concentration as much as possible. There is, however, a limit to the quantity of magnetic particles that can be added, as one can substitute a rare cell mixed in with other blood cells
20 for a rare cell mixed in with large quantities of magnetic particles upon separation. The latter condition does not markedly improve the ability to enumerate the cells of interest or to examine them.

[0036] The preferred magnetic particles for use in carrying out this invention are particles that behave as colloids. Such particles are characterized by their sub-micron particle size, which is generally less than about 200nm, and their stability to gravitational separation from solution for extended periods of time. In addition to the many other advantages,
25 this size range makes individual particles essentially invisible to analytical techniques commonly applied to cell analysis. Particles within the range of 90-150nm and having between 70-90% magnetic mass are contemplated for use in the present invention. Suitable magnetic particles are composed of a crystalline core of superparamagnetic material surrounded by molecules which are bonded, e.g., physically absorbed or covalently attached, to the magnetic core and which confer stabilizing colloidal properties. The coating material should preferably be applied in an amount effective to
30 prevent non-specific interactions between biological macromolecules found in the sample and the magnetic cores. Such biological macromolecules may include carbohydrates such as sialic acid residues on the surface of non-target cells, lectins, glycoproteins, and other membrane components. In addition, the material should contain as much magnetic mass per nanoparticle as possible. The size of the magnetic crystals comprising the core is sufficiently small that they do not contain a complete magnetic domain. The size of the nanoparticles is sufficiently small such that their Brownian energy
35 exceeds their magnetic moment. As a consequence, North Pole, South Pole alignment and subsequent mutual attraction/repulsion of these colloidal magnetic particles does not appear to occur even in moderately strong magnetic fields, contributing to their solution stability. Finally, the magnetic particles should be separable in high magnetic gradient external field separators. That characteristic facilitates sample handling and provides economic advantages over the more complicated internal gradient columns loaded with ferromagnetic beads or steel wool. Magnetic particles having
40 the above-described properties can be prepared by modification of base materials described in U.S. Patents #4,795,698, #5,597,531, and #5,698,271.

[0037] An improved method for making particles is described in U.S. Patent #5,698,271. These materials are an improvement over those disclosed in the '531 patent in that the process includes a high temperature coating step which markedly increases the level of coating. Nanoparticles made with bovine serum albumin (BSA) coating using this process,
45 for example, have a 3-5-fold lower non-specific binding characteristic for cells when compared to the DC-BSA materials of '531. This decrease in non-specific binding has been shown to be directly due to the increased level of BSA coating material. When such nanoparticles were treated so as to remove BSA coating, non-specific binding returns to high levels. It was thus determined that a direct relationship exists between the amount of BSA coated on iron oxide crystal surfaces and the nonspecific binding of cells. Typically, the non-specific binding of cells from whole blood with these particles
50 was 0.3%, which is significantly better than those, produced from '531. Thus, from 10ml of whole blood there would be about 200,000 non-target cells that would also be isolated with the cells targeted for enrichment.

[0038] Since small nanoparticles (30-70nm) will diffuse more readily they will preferentially label cells compared with their larger counterparts. When very high gradients are used, such as in internal gradient columns, the performance of these materials, regardless of size, makes little difference. On the other hand, when using external gradients, or gradients of lesser magnitude than can be generated on microbead or steel wool columns, the occupancy of small nanoparticles
55 on cells has a significant effect. This was conclusively shown to be the case by fractionating DC nanoparticles and studying the effects on recovery. Based on these studies and other optimization experiments, means for fractionating nanoparticles magnetically or on columns was established where base coated magnetic particles could be prepared

EP 1 425 294 B1

that were devoid of excessively small or large nanoparticles. For example, base coated particles of mean diameter 100nm can be produced which contain at best trace amounts of material smaller than 80nm or over 130nm. Similarly material of about 120nm can be made with no appreciable material smaller than 90-95nm and over 160nm. Such materials performed optimally with regard to recovery and could be made sub-optimal by the inclusion of 60-70nm nanoparticles. The preferred particle size range for use in practicing this invention is 90-150nm for base coated magnetic particles, e.g., BSA-coated magnetite.

[0039] Based on the foregoing, high gradient magnetic separation with an external field device employing highly magnetic, low non-specific binding, colloidal magnetic particles is the method of choice for separating a cell subset of interest from a mixed population of eukaryotic cells, particularly if the subset of interest comprises but a small fraction of the entire population. Such materials, because of their diffusive properties, readily find and magnetically label rare events, such as tumor cells in blood. For magnetic separations for tumor cell analysis to be successful, the magnetic particles must be specific for epitopes that are not present on hematopoietic cells.

[0040] A large variety of analytical methods and criteria are used to identify tumor cells, and the first attempts are being undertaken to standardize criteria that define what constitutes a tumor cell by immunocytochemistry. In this study, blood samples from prostate cancer patients were immunomagnetically enriched for cells that expressed EpCAM. Tumor cells were identified by the expression of the cytoskeletal proteins cytokeratin (CK+), the absence of the common leukocyte antigen CD45 (CD45-) and the presence of nucleic acids (NA+) by multicolor fluorescence analysis. Rare events or rare cells can be immunophenotyped by both flowcytometry and fluorescence microscopy. Flowcytometric analysis excels in its ability to reproducibly quantify even low levels of fluorescence whereas microscopy has the better specificity as morphological features can aid in the classification of the immunophenotypically identified objects. Although there was a correlation between the number of CTC detected in blood of prostate cancer patients by flowcytometry and microscopy, microscopic examination of the CK+, CD45-, NA+ objects showed that only few of the objects appeared as intact cells. This observation agrees with other reports that showed apoptosis in a substantial portion of circulating tumor cells.

[0041] The terms "biological specimen" or "biological sample" may be used interchangeably, and refer to a small portion of fluid or tissue taken from a human test subject that is suspected to contain cells of interest, and is to be analyzed. A biological specimen refers to the fluidic portion, the cellular portion, and the portion containing soluble material. Biological specimens or biological samples include, without limit bodily fluids, such as peripheral blood, tissue homogenates, nipple aspirates, colonic lavage, sputum, bronchial (alveolar) lavage, pleural fluids, peritoneal fluids, pericardial fluids, urine, and any other source of cells that is obtainable from a human test subject. An exemplary tissue homogenate may be obtained from the sentinel node in a breast cancer patient.

[0042] The term "rare cells" is defined herein as cells that are not normally present in biological specimens, but may be present as an indicator of an abnormal condition, such as infectious disease, chronic disease, injury, or pregnancy. Rare cells also refer to cells that may be normally present in biological specimens, but are present with a frequency several orders of magnitude less than cells typically present in a normal biological specimen.

[0043] The term "determinant", when used in reference to any of the foregoing target bioentities, refers broadly to chemical mosaics present on macromolecular antigens that often induce an immune response. Determinants may also be used interchangeably with "epitopes". A "biospecific ligand" or a "biospecific reagent," used interchangeably herein, may specifically bind determinants. A determinant refers to that portion of the target bioentity involved in, and responsible for, selective binding to a specific binding substance (such as a ligand or reagent), the presence of which is required for selective binding to occur. In fundamental terms, determinants are molecular contact regions on target bioentities that are recognized by agents, ligands and/or reagents having binding affinity therefor, in specific binding pair reactions.

[0044] The term "specific binding pair" as used herein includes antigen-antibody, receptor-hormone, receptor-ligand, agonist-antagonist, lectin-carbohydrate, nucleic acid (RNA or DNA) hybridizing sequences, Fc receptor or mouse IgG-protein A, avidin-biotin, streptavidin-biotin and virus-receptor interactions.

[0045] The term "detectable label" is used to herein to refer to any substance whose detection or measurement, either directly or indirectly, by physical or chemical means, is indicative of the presence of the target bioentity in the test sample. Representative examples of useful detectable labels, include, but are not limited to the following: molecules or ions directly or indirectly detectable based on light absorbance, fluorescence, reflectance, light scatter, phosphorescence, or luminescence properties; molecules or ions detectable by their radioactive properties; molecules or ions detectable by their nuclear magnetic resonance or paramagnetic properties. Included among the group of molecules indirectly detectable based on light absorbance or fluorescence, for example, are various enzymes which cause appropriate substrates to convert (e.g., from non-light absorbing to light absorbing molecules, or from non-fluorescent to fluorescent molecules). Analysis can be performed using any of a number of commonly used platforms, including multiparameter flow cytometry, immunofluorescent microscopy, laser scanning cytometry, bright field base image analysis, capillary volumetry, spectral imaging analysis, manual cell analysis, CellSpotter® analysis, CellTracks™ analysis, and automated cell analysis.

[0046] The phrase "to the substantial exclusion of" refers to the specificity of the binding reaction between the biospecific

EP 1 425 294 B1

ligand or biospecific reagent and its corresponding target determinant. Biospecific ligands and reagents have specific binding activity for their target determinant yet may also exhibit a low level of non-specific binding to other sample components.

[0047] The phrase "early stage cancer" is used interchangeably herein with "Stage I" or "Stage II" cancer and refers to those cancers that have been clinically determined to be organ-confined. Also included are tumors too small to be detected by conventional methods such as mammography for breast cancer patients, or X-rays for lung cancer patients. While mammography can detect tumors having approximately 2×10^8 cells, the methods of the present invention should enable detection of circulating cancer cells from tumors approximating this size or smaller.

[0048] The term "enrichment" as used herein refers to the process of substantially increasing the ratio of target bioentities (e.g., tumor cells) to non-target materials in the processed analytical sample compared to the ratio in the original biological sample. In cases where peripheral blood is used as the starting materials, red cells are not counted when assessing the extent of enrichment. Using the method of the present invention, circulating epithelial cells may be enriched relative to leucocytes to the extent of at least 2,500 fold, more preferably 5,000 fold and most preferably 10,000 fold.

[0049] The terms "anti-coagulant" or "anti-coagulating agent" may be used interchangeably, and refer to compositions that are added to biological specimens for the purpose of inhibiting any undesired natural or artificial coagulation. An example of coagulation is blood clotting and common anti-coagulants are chelating agents, exemplified by ethylenediamine tetraacetic acid (EDTA), diethylenetriamine pentaacetic acid (DTPA), 1,2-diaminocyclohexane tetraacetic acid (DCTA), ethylenebis(oxyethylenenitrilo) tetraacetic acid (EGTA), or by complexing agents, such as heparin, and heparin species, such as heparin sulfate and low-molecular weight heparins. This may be further collectively defined as "clumping" or "clump formation". However, such clumps must be differentiated from "clusters" or aggregates of CTC that are counted as a single Intact CTC if they meet the classification criteria for Intact CTC.

[0050] Clusters of CTC are believed to have greater proliferative potential than single CTC and their presence is thus diagnostically highly significant. One possible cause for an increased propensity to establish secondary metastatic tumor sites may be the virtue of their adhesiveness. An even more likely cause is the actual size of a CTC cluster; larger clusters will become lodged in small diameter capillaries or pores in bone. Once there, the viability of the cells in the cluster would determine the chance of survivability at the new metastatic site.

[0051] The ideal "stabilizer" or "preservative" (herein used interchangeably) is defined as a composition capable of preserving target cells of interest present in a biological specimen, while minimizing the formation of interfering aggregates and cellular debris in the biological specimen, which in any way can impede the isolation, detection, and enumeration of targets cells, and their differentiation from non-target cells. In other words, when combined with an anti-coagulating agent, a stabilizing agent should not counteract the anti-coagulating agent's performance. Conversely, the anti-coagulating agent should not interfere with the performance of the stabilizing agent. Additionally, the disclosed stabilizers also serve a third function of fixing, and thereby stabilizing, permeabilized cells, wherein the expressions "permeabilized" or "permeabilization" and "fixing", "fixed" or "fixation" are used as conventionally defined in cell biology. The description of stabilizing agents herein implies using these agents at appropriate concentrations or amounts, which would be readily apparent to one skilled in cell biology, where the concentration or amount is effective to stabilize the target cells without causing damage. One using the compositions, methods, and apparatus of this invention for the purpose of preserving rare cells would obviously not use them in ways to damage or destroy these same rare cells, and would therefore inherently select appropriate concentrations or amounts. For example, the formaldehyde donor imidazolidinyl urea has been found to be effective at a preferred concentration of 0.1-10%, more preferably at 0.5-5% and most preferably at about 1-3% of the volume of said specimen. An additional agent, such as polyethylene glycol has also been found to be effective, when added at a preferred concentration of about 0.1 % to about 5%, more preferably about 0.1 % to about 1%, and most preferably about 0.1% to about 0.5% of the specimen volume.

[0052] A stabilizing agent must be capable of preserving a sample for at least a few hours. However, it has been shown that samples can be stabilized for at least up to 72 hours. Such long-term stability is important in cases where the sample is obtained in a location that is distant to the location where processing and analysis will occur. Furthermore, the sample must be stabilized against mechanical damage during transport.

[0053] Stabilizing agents are necessary to discriminate between *in vivo* tumor cell disintegration and disintegration due to *in vitro* sample degradation. Therefore, stabilizing agent compositions, as well as methods and apparatus for their use, are described in WO 03/018757, published on 6th March 2003.

[0054] The terms "Obvious cells" or "intact cells" may be used interchangeably, and refer to cells found during imaging analysis that contain nucleic acid and cytokeratin. These cells are usually visually round or oval, but may sometimes be polygonal or elongated, and appear as individual cells or clusters of cells. The nucleic acid area (i.e. labeled by nucleic acid dye) is smaller than the cytoplasmic area (i.e. labeled by anti-cytokeratin), and is surrounded by the cytoplasmic area.

[0055] The terms "suspicious cells", "Suspect cells", or "fragments" may be used interchangeably, and refer to cells found during imaging analysis that resemble intact cells, but are not as visually distinct as intact cells. Based on imaging analysis, there are a number of possible types of Suspect cells, including:

EP 1 425 294 B1

1. Enucleated cells, which are shaped like Obvious cells, are positively stained for cytokeratin, but negative for nucleic acid;
2. Speckled or punctate cells, which are positively stained for nucleic acid, but have irregularly-stained cytokeratin; and
3. Amorphic cells, which stain positively for cytokeratin and nucleic acid, but are irregular in shape, or unusually large.

These suspicious cells are of interest in this invention because they may give additional information to the nature of the CTC, as well as the patient's disease. It is possible that staining or image artifacts may be observed during analysis. For example, enucleated cells sometimes appear to have a "ghost" region where the nucleus should have stained, but the corresponding region is nucleic acid negative. This may be caused by a number of external factors, including the labeling or imaging techniques. Also, cells have been observed with "detached" nuclei. While this may possibly indicate a cell releasing its nucleus, it is more likely that this appears due to an artifact of the imaging system. However, such "artifacts," when real, give valuable information about what may be happening to the intact cells. Therefore, as part of this invention, suspicious cells will be more closely analyzed.

[0056] Cell fragments are different than "debris" in that debris does not necessarily resemble a cell. The term debris as used herein, refers to unclassified objects that are specifically or non-specifically labeled during processing, and are visible as images during analysis, but are distinct from intact and/or Suspect cells. For example, it has been observed that damaged cells will release nuclear material. During processing, this nuclear material may be non-specifically magnetically labeled, and subsequently labeled with the nucleic acid stain. During analysis, the magnetically labeled and stained nuclear material can be observed when it has cytokeratin still attached. There are other objects that are similarly magnetically selected and stained which appear during analysis that are classified as debris.

[0057] The term "morphological analysis" as used herein, refers to visually observable characteristics for an object, such as size, shape, or the presence/absence of certain features. In order to visualize morphological features, an object is typically non-specifically stained. The term "epitope analysis" as used herein, refers to observations made on objects that have been labeled for certain epitopes. In order to visualize epitopic features, an object is typically specifically stained or labeled. Morphological analysis may be combined with epitope analysis to provide a more complete analysis of an object.

[0058] The importance of further visual observation is apparent when fragments and debris are often classified as "Not Assigned Events," or "Unassigned events." These terms arise from non-visual analysis, such as with flow cytometry. Because flow cytometry does not image objects, any event not falling in the specified populations that meet the criteria for the target cells, or the non-target cells (as is the case when non-specifically carried over WBC are negatively labeled), will fall outside either of these populations. However, as will be apparent throughout this specification, these unassigned events are important.

[0059] Figure 1 is a model of various CTC stages, including shedding and metastasis. Figure 1a. shows these stages for cells, clusters, fragments, and debris. Figure 1b. shows actual images from samples at these same stages. The images of cells clusters, fragments, and debris were taken from CellSpotter® analyses of patient samples. The images of tissue samples (Origin and Metastatic sites) were taken from elsewhere (Manual of Cytology, American Society of Clinical Pathologists Press. 1983).

[0060] Briefly, a single cell shed from a primary tumor into the blood either survives or dies in blood. If it survives, it may possibly divide in blood, or colonize at a secondary site. If the cell dies, depending on the method, the cell degrades into various types of fragments or debris. Another possibility is a cluster of cells is shed from a primary tumor into the blood, where it may dissociate into single cells, or remain intact, and colonize at a secondary site. If the cluster dissociates, it can behave similar to the single cell described above. If the cluster remains intact, it is more likely to for a secondary colony for the reasons described above, which includes the large diameter cluster becoming lodged in a small diameter capillary. Once lodged, if the cells are viable, the cluster would form a new tumor.

[0061] The presence of fragments and debris with very few intact cells suggests that there will be little chance of metastasis. Fragmented cells will not divide, and cannot form secondary tumors. Indeed, only intact CTC or possibly CTC clusters would be capable of colonizing secondary sites. Identification of antigens that play a role in the adhesion and penetration process may help. Follow up and assessment of metastatic sites of the patients with and without clusters will also provide further insight. Nuclear morphology is used to determine the activity status and abnormality of a cell. Chromatin clumping, the presence or absence of nucleoli, and hyperchromasia, are criteria used to determine whether a cell is benign or malignant, reacting to a immune response, or reacting to treatment. The cytoplasmic morphology is used to determine the level of differentiation (i.e. tissue of origin). For example, cytoplasmic morphology can classify cells as squamous versus glandular.

[0062] During blood draw and subsequent specimen processing, the surviving battered tumor cells present in the peripheral circulation may be further stressed and damaged by turbulence during blood draw into an evacuated tube and by specimen processing, e.g. transport of the blood tube and mixing prior to analysis. Such mechanical damage is additional to on-going immunological, apoptotic, and necrotic processes leading to destruction of CTC that occur *in vitro*

EP 1 425 294 B1

in a time dependent manner. We have found that the longer the specimen is stored, the greater the loss of CTC, and the larger the amounts of interfering debris and/or aggregates. Indeed, data presented in this specification (Figures 2 and 3) show dramatic declines in CTC counts in several blood specimens stored at room temperature for 24hrs or longer, indicating substantial *in vitro* destruction of CTC after blood draw. While the losses of hematopoietic cells are well known phenomena and the subject of above-cited patents by Streck Labs and by others, the occurrence of mechanical damage due to mixing or transport have to date not been recognized factors in the loss of CTC. The formation of cellular debris and the interfering effects of accumulating debris and/or aggregates in the analysis of CTC have similarly been unrecognized to date. It appears to be most evident and problematic in highly sensitive enrichment assays requiring processing of relatively large blood volumes (5-50mL), and subsequent microscopic detection or imaging of target cells after volume reduction (less than 1mL). Such debris are either not normally seen, or do not interfere in conventional non-enrichment assays, for example, by flow cytometry or in enrichment by density gradients methods.

[0063] To explore if these damaged epithelial cells and epithelial cell fragments observed in patients could be caused by apoptosis of tumor cells induced by chemotherapy, a model to mimic tumor cell death was developed. Cells of the prostate cell line LnCaP were cultured with or without paclitaxel and spiked into blood of healthy donors. The immunomagnetically selected cells of the paclitaxel treated samples analyzed by CellSpotter® resembled those observed in the patient blood samples. Cells treated with paclitaxel displayed signs of apoptosis. The punctate cytokeratin staining pattern of the cells appear to correspond with a collapse of the cytoskeletal proteins (Figure 4B vs. 6B). The initiating event in the sequence resulting from the microtubule stabilizing effects of paclitaxel which in turn may activate the pro-apoptotic gene Bim that senses cytoskeletal distress. Further evidence of caspase-cleaved cytokeratin resulting from apoptosis was obtained with the M30 Cytodeath antibody (Roche Applied Science, Mannheim, Germany) that recognizes an epitope of cytokeratin 18 that is only exposed following caspase cleavage in early apoptosis. Only the paclitaxel treated LnCaP cells stained with M30 and most of the dimmer cytokeratin cells stained with M30, which is consistent with cells undergoing apoptosis.

[0064] It should be noted that a number of different cell analysis platforms can be used to identify and enumerate cells in the enriched samples. Examples of such analytical platforms are Immunicon's CellSpotter® system, a magnetic cell immobilization and analysis system, using microscopic detection for manual observation of cells described in Example II, and the CellTracks™ system, an a more advanced automatic optical scanning system. These two analytical platforms are described in US Patents #5,876,593; #5,985,153 and #6,136,182, as disclosing the respective apparatus and methods for manual or automated quantitative and qualitative cell analysis.

[0065] Other analysis platforms include laser scanning Cytometry (Compucyte), bright field base image analysis (Chromavision), and capillary Volumetry (Biometric Imaging).

[0066] The enumeration of circulating epithelial cells in blood using the methods and compositions of a preferred embodiment of the present invention is achieved by immunomagnetic selection (enrichment) of epithelial cells from blood followed by the analysis of the samples. The immunomagnetic sample preparation is important for reducing sample volume and obtaining as much as a 10^4 fold enrichment of the target (epithelial) cells. The reagents used for the multi-parameter flow cytometric analysis are optimized such that epithelial cells are located in a unique position in the multi-dimensional space created by the listmode acquisition of two light scatter and three fluorescence parameters. These include

1. an antibody against the pan-leukocyte antigen, CD45 to identify leucocytes (non-tumor cells);
2. a cell type specific or nucleic acid dye which allows exclusion of residual red blood cells, platelets and other non-nucleated events; and
3. a biospecific reagent or antibody directed against cytokeratin or an antibody having specificity for an EpCAM epitope which differs from that used to immunomagnetically select the cells.

[0067] It will be recognized by those skilled in the art that the method of analysis of the enriched tumor cell population will depend on the intended use of the invention. For example, in screening for cancers or monitoring for recurrence of disease, as described hereinbelow, the numbers of circulating epithelial cells can be very low. Since there is some "normal" level of epithelial cells, (very likely introduced during venipuncture), a method of analysis that identifies epithelial cells as normal or tumor cells is desirable. In that case, microscopy based analyses may prove to be the most accurate. Such examination might also include examination of morphology, identification of known tumor diathesis associated molecules (e.g., oncogenes).

Patients

[0068] Patients' age range was 47-91 year (mean 74), with initial diagnosis 2 to 10 years prior to study. Medical records were reviewed for therapy and stage. Patients and healthy volunteers signed an informed consent under an approved research study. Blood was drawn into 10ml EDTA Vacutainer™ tubes (Becton-Dickinson, NJ). Samples were kept at

EP 1 425 294 B1

room temperature and processed within 6 hours after collection unless indicated otherwise.

Sample Preparation

[0069] Magnetic nanoparticles labeled with monoclonal antibodies identifying epithelial cell adhesion molecule (Ep-CAM) were used to label and separate by magnetic means epithelial cells from hematopoietic cells as taught in commonly-owned US Patent #6,365,362, and US 2002/0172987 A1, published on 21st November 2002. The magnetically captured cells resuspended in a volume of 200 μ l are fluorescently labeled to differentiate between hematopoietic and epithelial cells. A monoclonal antibody that recognizes keratins 4, 5, 6, 8, 10, 13, and 18, conjugated to Phycoerythrin (CK-PE) was used to identify epithelial cells and a monoclonal antibody that recognizes CD45 was used to identify leukocytes and identify hematopoietic cells that non-specifically bind to cytokeratin.

[0070] For multicolor fluorescent microscopy (CellSpotter[®]) analysis CD45 was conjugated to Allophycocyanin (CD45-APC, Caltag, CA) whereas for flow cytometric analysis peridinin chlorophyll protein conjugated CD45 (CD45-PerCP, BDIS San Jose, CA) was used. The nucleic acid specific dye DAPI (4,6-diamidino-2-phenylindole) was used to identify and visualize the nucleus in the CellSpotter[®] system and the nucleic acid dye in the Procount system (BDIS, San Jose, CA) was used to identify cells by flow cytometry. After incubation, the excess staining reagents were aspirated and the captured cells were resuspended and transferred into a 12x75 mm tube for flow cytometric analysis or to a CellSpotter[®] analysis chamber (as described in US 2002/0109838 A1, published on 15th August 2002) contained within a magnetic yoke assembly that holds the chamber between two magnets (Captivate, Molecular Probes, OR).

Example 1**Sample Analysis via Flow Cytometry**

[0071] Samples were analyzed on a FACSCalibur flow cytometer equipped with a 488nm Argon ion laser (BDIS, San Jose, CA). Data acquisition was performed with CellQuest (BDIS, San Jose, CA) using a threshold on the fluorescence of the nucleic acid dye. The acquisition was halted after 8000 beads or 80% of the sample was analyzed. Multiparameter data analysis was performed on the listmode data (Paint-A-Gate^{Pro}, BDIS, San Jose, CA). Analysis criteria for CTC events included size defined by forward light scatter, granularity defined by orthogonal light scatter, positive staining with the PE-labeled anti-cytokeratin MAb and no staining with the PerCP-labeled anti-CD45 MAb. For each sample, the number of events present in the region typical for epithelial cells was multiplied by 1.25 to account for the sample volume not analyzed by flow cytometry.

[0072] Figure 2 Panels A, B and C shows the flow cytometric analysis of a blood sample of a patient with metastatic prostate cancer. Two vertical lines in Panel B illustrate the low and high boundary of nucleic acid (NAD) content of leukocytes (red dots). CTC candidates express Cytokeratin (CK+), lack CD45 (CD45-) and contain nucleic acids (NAD+). CTC candidates having NAD equal or higher than leukocytes are considered cells and are depicted black. CK+, CD45- events with NAD content less than leukocytes were not considered target cells and depicted blue. The blue events were clearly smaller as compared with the black colored CTC as evident by the smaller forward light scatter signals. The threshold on the NAD staining intensity clearly excluded a large portion of CK+, CD45- events with even lower NAD staining intensity. In analysis of blood samples from healthy donors few such CK+, CD45- events are observed suggesting that this phenomenon is related to cancer. A typical example of an analysis of a blood from a healthy donor is shown in Figures 2D, 2E, and 2F.

Example 2**Sample Analysis via CellSpotter[®]**

[0073] The CellSpotter[®] system consists of a microscope with a Mercury Arc Lamp, a 10X objective, a high resolution X, Y, Z stage and a four-filter cube changer. Excitation, dichroic and emission filters in each of four cubes were for DAPI 365nm/400nm/400nm, for DiOC16 480nm/ 495nm/ 510nm, for PE 546nm/ 560nm/ 580nm and for APC 620nm/ 660nm/ 700nm. Images were acquired with a digital camera connected to a digital frame grabber. The surface of the chamber is 80.2 mm² and 4 rows of 35 images for each of the 4 filters resulting in 560 images have to be acquired to cover the complete surface. The CellSpotter[®] acquisition program automatically determines the region over which the images are to be acquired, the number of images to acquire, the position of each image and the microscope focus to use at each position. All the images from a sample are logged into a directory that is unique to the specific sample identification. An algorithm is applied on all of the images acquired from a sample to search for locations that stain for DAPI and CK-PE. If the staining area is consistent with that of a potential tumor cell (DAPI+, CK-PE+) the software stores the location of these areas in a database. The software displays thumbnails of each of the boxes and the user can confirm that the

EP 1 425 294 B1

images represented in the row are consistent with tumor cells, or stain with the leukocyte marker CD45. The software tabulates the checked boxes for each sample and the information is stored in the database.

[0074] Figure 3 shows examples of CellSpotter® analysis of a blood sample from a patient with metastatic prostate cancer. Regions that potentially contain tumor cells are displayed in rows of thumbnails. The ruler in the left lower corner of the figure indicates the sizes of the thumbnails. From right to left these thumbnails represent nuclear (DAPI), cytoplasmic cytokeratin (CK-PE), control cells stained with a membrane dye (DiOC₁₆(3)) and surface CD45 (CD45-APC) staining. The composite images shown at the left show a false color overlay of the purple nuclear (DAPI) and green cytoplasmic (CK-PE) staining. The check box beside the composite image allow the user to confirm that the images represented in the row are consistent with tumor cells and the check box beside the CD45-APC image is to confirm that a leukocyte or tumor cell stain non-specifically. In this patient sample, the software detected 2761 rows of thumbnails that demonstrated staining consistent with tumor cells. Eighteen of the 2761 rows are shown in the figure labeled 1631-1640 and 1869-1876. Rows numbered 1631, 1636, 1638, 1640, and 1873-1876 are checked off and display features of CTC defined as a size greater than 4µm, the presence of a nucleus surrounded by cytoplasmic cytokeratin staining and absence of DiOC₁₆(3) and CD45 staining. Note the difference in appearance of the tumor cells: the cell in row 1638 is large and the one in row 1640 is significantly smaller. The immunophenotype of the events in rows 1634 and 1869 are consistent with tumor cells but their morphology is not consistent with intact cells. The thumbnails in row 1869 shows a large nucleus and speckled cytoplasmic due to retraction of cytoskeletal proteins consistent with apoptosis of the cell. The thumbnail in row 1634 shows a damaged cell that appears to extrude its nucleus. The thumbnail shown in row 1632 shows a cell that stains both with cytokeratin as well as CD45 and is either a tumor cell non-specifically binding to CD45 or a leukocyte non specifically staining with cytokeratin. In this instance the morphology of the cell closely resembles that of a lymphocyte. The thumbnails shown in rows 1633, 1635, 1637, 1639, 1870 and 1872 shows cytokeratin staining objects that are larger than 4 µm but have no resemblance to cells. The cytokeratin staining objects in thumbnails 1637, 1639 and 1872 are in close proximity of a leukocyte.

[0075] Based on observation of images of CTC candidates in several patient samples, CTC were classified into three categories: intact CTC, damaged CTC, and CTC fragments all not staining with CD45 and not appearing in the DiOC₁₆(3) filter. Figure 4 displays examples of the three categories of CTC isolated from a single tube of blood of a patient with metastatic prostate cancer undergoing therapy. Intact tumor cells shown in Figure 3A were defined as objects larger than 4µm with a relatively smooth cytoplasmic membrane, cytoskeletal proteins throughout the cytoplasm, and an intact nucleus encompassed within the nucleus. Damaged CTC shown in Figure 4B were defined as objects larger than 4µm with speckled cytokeratin staining or ragged cytoplasmic membrane, and a nucleus associated with the cytokeratin staining. Tumor cell fragments shown in Figure 4C were defined as round cytokeratin staining objects larger than 4µm with or without association of nuclear material that had no morphological resemblance to a cell.

Example 3

CTC in Prostate Cancer Patients

[0076] CTC were enumerated in 18 blood samples of prostate cancer patients and 27 samples from healthy individuals by both flow cytometry and CellSpotter®. The results shown in Table 1 were sorted by increasing number of intact CTC detected by CellSpotter®.

Table 1- Enumeration of CTC by CellSpotter® and flow cytometry in 18 blood samples of prostate cancer patients and 27 samples from healthy individuals.

Patient Sample	CellSpotter®						Flow Cytometry
	Intact CTC		Suspect CTC		Not Assigned Events		CK+CD45-NA+
	#	%	#	%	#	%	#
1	0	0	1	50	1	50	5
2	0	0	2	100	0	0	12
3	0	0	2	66	1	34	1
4	0	0	2	50	2	50	0
5	0	0	2	29	5	71	5

EP 1 425 294 B1

(continued)

Patient Sample	CellSpotter®						Flow Cytometry
	Intact CTC		Suspect CTC		Not Assigned Events		CK+CD45-NA+
	#	%	#	%	#	%	#
6	0	0	3	60	2	40	18
7	0	0	3	38	5	62	0
8	0	0	7	44	9	56	10
9	0	0	13	76	4	24	2
10	1	5		5	20	90	4
11	1	10	4	40	5	50	0
12	2	22	1	11	6	67	4
13	28	6	7	1	441	93	69
14	70	5	168	12	1204	83	683
15	322	3	448	13	4244	87	500
16	350	5	112	2	5924	93	723
17	350	2	1429	9	14412	89	2420
18	742	17	112	2	3641	81	310
Mean	-	4%	-	34%	-	62%	-
27 samples from healthy donors							
Mean	0.04		0.96		4.96		0.7
SD	0.19		1.85		3.98		1.14
Min	0		0		0		0
Max	1		7		15		4
# - number CTC in 7.5 ml blood % - percentage of all CTC detected by CellSpotter®							

[0077] In the CellSpotter® analysis, the proportion of intact CTC clearly constituted the smallest fraction of CTC and ranged from 0% to 22% of all CTC (mean 4%). The proportion of damaged CTC ranged from 1 % to 100% (mean 34%) and the CTC fragments constituted the largest portion of CTC ranging from 0% to 93% (mean 62%). The distribution of CTC over the three categories between the patients varied considerably as amplified by a lack of correlation between intact CTC and damaged CTC ($R^2 = 0.20$) and intact CTC and CTC fragments ($R^2 = 0.42$) and some correlation between damaged CTC and CTC fragments ($R^2 = 0.88$). Comparison of intact CTC by CellSpotter® and CTC enumerated by flow cytometry showed no significant correlation ($R^2 = 0.26$) whereas significant correlations were found between the damaged CTC and CTC by flow cytometry ($R^2 = 0.92$) and CTC fragments and CTC by flow cytometry ($R^2 = 0.93$). Comparison of the CTC detected by flow cytometry and CellSpotter® suggests that CTC detected by flow cytometry encompass intact CTC as well as damaged CTC and to a certain extent, CTC fragments.

Example 4

Mimicking cell damage by in-vitro induction of apoptosis in LnCaP cells

[0078] To investigate the effect of apoptosis induced by cytotoxic agents on flow cytometric and CellSpotter® analysis of CTC, cells from the prostate cell line LnCaP were cultured in the presence or absence of 40nM paclitaxel for 72 hours. Following incubation, untreated LnCaP cells demonstrated a viability of >95% by trypan blue exclusion and 33% for the paclitaxel treated cells. The treated and untreated LnCaP cells were spiked into blood of healthy donors, selected by the ferrofluid methods described above, and analyzed by the CellSpotter® system. In experiments in which LnCaP cells were spiked into blood that were not treated with paclitaxel greater than 95% of the LnCaP cells were classified as intact tumor cells. The morphologic appearance of the paclitaxel treated LnCaP cells showed close resemblance to that of the

EP 1 425 294 B1

CTC observed in the patient samples and are shown in Figure 6. Intact LnCaP cells that survived paclitaxel treatment are shown in Figure 6A, damaged LnCaP, of which the majority show speckled cytokeratin staining, are shown in Figure 6B, and tumor fragments are shown in Figure 6C.

[0079] Normal blood samples spiked with paclitaxel treated and untreated LnCaP cells were also prepared for flow cytometric analysis. In Figures 2G, 2H, and 2I, the flow cytometric analysis of a blood sample spiked with 501 LnCaP cells is shown. A predominantly bright cytokeratin positive population with a nucleic acid content greater than normal human leukocytes and relatively large size as illustrated by the large forward light scatter signals are shown and depicted black in the figure. Only few CK+, CD45- events with NAD content less than leukocytes and depicted blue are detected in the sample. Figures 2J, 2K, and 2L shows the flow cytometric analysis of paclitaxel treated LnCaP cells spiked in blood. In contrast to viable LnCaP cells, a wide distribution of cytokeratin staining was observed with a significant portion of the population demonstrating a decreased concentration of nucleic acid content. In addition, numerous small cytokeratin positive events with less nucleic acid content as leukocytes were observed. The pattern of the patient closely resembled that of the pattern of the paclitaxel treated LnCaP cells supporting the hypothesis that the CTC detected by flow cytometry represent intact CTC as well as a variety of disintegrating cells in blood of cancer patients.

[0080] The data shown above demonstrate that in the blood of patients with prostate cancer, CTC detected by both flow cytometry and CellSpotter® are comprised of intact cells and cells of cells at various stages of disintegration. The apoptosis induced *in vitro* by paclitaxel suggests that the detected CTC in patient blood samples are undergoing apoptosis, necrosis, or *in vivo* damage to a varying degree caused by the treatment or therapy, mechanical damage by passage through the vascular system, or by the immune system.

[0081] Another source of cell disintegration, caused *in vitro* could however, be introduced by the sample preparation or the lack of stabilization of CTC or other blood components after blood draw. To investigate the effect of sample aging, known to cause damage, blood samples drawn from 12 patients with prostate cancer were processed and analyzed by flow cytometry within two hours, after 24 hours, and after 6 and 18 hours if sufficient blood was available. In 8 of the 12 patient samples, CTC were detected at a level greater than the mean +3SD of that detected in normal donors. As shown in Table 2, a loss of CTC with sample aging was observed in all 8 samples.

Table 2 - Enumeration of CTC by flow cytometry in 8 blood samples of prostate cancer patients processed and analyzed at different time points after blood draw

Patient #	Time after blood draw			
	< 2hr #CTC	~6hr #CTC	~18hr #CTC	~24hr #CTC
1	5	-	-	0
2 8 9 2				3
3	15	-	-	0
4	31	-	-	3
5	44	-	-	8
6	45	-	-	1
7	49	38	19	26
8	78	-	-	0
hr = hours #CTC = number of CTC in 5 ml blood				

[0082] Significant reductions in the number of CTC were detected when blood processing was delayed demonstrating the fragility of CTC, and making it necessary to process non-stabilized blood samples no later than six hours after blood draw to obtain accurate CTC counts. To reliably assess if clinically relevant information is contained within the different stages of tumor cell degradation, a blood preservative is needed that stabilizes CTC at the time of blood draw to obtain an accurate reflection of what is occurring inside the body. Furthermore, the sample preparation method for sensitive assays used to enrich for CTC requires that all classes of CTC are captured, and therefore excludes the use of traditional density gradient separation methods in the prior art.

EP 1 425 294 B1

Example 5

Obvious CTC and Suspect CTC are important indicators

5 [0083] It is important to be able to distinguish between *in vivo* and *in vitro* damage for sensitive assays, such as those described here. This is especially evident when the assay attempts to determine the effectiveness of treatments or therapies, which are known to cause *in vivo* cellular damage. If sample handling, processing, or analysis were to result in damaging the target cells, forming Suspect cells, fragments, or debris, the assay will not give meaningful results.

10 [0084] An assay was used to directly detect CTC in 100 μ l of blood without any enrichment method by flow cytometry. The 100 μ l assay detects only EpCAM positive cells and the sensitivity is very low. However, some advanced stage cancer patients with high CTC counts are expected to be observable. This assay should give a reliable confirmatory estimation of CTC because it is a direct assay that involves no manipulation. Data were generated with several patient samples using the assay to answer several questions.

15 [0085] The 100 μ l assay categorizes cells based on properties such as size and staining intensity. Obvious CTC have bright nucleic acid staining (similar to leukocytes), positive EpCAM antigen staining and size similar to leukocytes or larger. Suspect CTC are any objects positive for EpCAM antibody but not characterized as Obvious CTC (i.e. dim nucleic acid, size smaller than leukocytes). The assay identifies objects from both categories.

20 [0086] Figure 5 shows the presence of Obvious and Suspect CTC in blood as determined by the 100 μ l assay. The Suspect CTC are not created during sample processing (*in vitro* damage) as the 100 μ l assay is a direct assay and does not involve any separation or wash steps. The data above also show there is a relationship between the number of Obvious and Suspect CTC. The number of Suspect CTC seems to increase as the number of Obvious CTC increases. When the numbers of Suspect versus Obvious CTC is plotted, the slope of 2.92 indicates the proportion of Suspect CTC present in sample when compared to Obvious CTC. The correlation coefficient of $r^2 = 0.97$ shows an excellent correlation between Obvious CTC and Suspect CTC for a number of clinical samples. In addition, Suspect CTC are also seen in the ferrofluid-selection assay, and have properties similar to Suspect CTC detected in the blood by the direct assay. It is important to include Suspect CTC in addition to Obvious CTC in total tumor cell count.

25 [0087] An important question is how the data from the 100 μ l assay compares with ferrofluid-selected CTC (enriched CTC). Does the ferrofluid assay quantitatively detect CTC? Another question is what is the recovery of CTC in the ferrofluid-selected assay if the flow assay data is correct. The three main factors determining the recovery of CTC in the 100 μ l assay are:

- EpCAM density,
- cytokeratin positivity, and
- nucleus positivity.

30 The Suspect CTC have lower EpCAM density compared to Obvious CTC and significance of this is not yet well understood.

35 [0088] A comparison was made of Obvious and Suspect CTC by the 100 μ l assay to the ferrofluid-selection assay using 7ml of blood. This data was obtained from prostate patient samples and analyzed by flow cytometry. Both Obvious and Suspect CTC increased with storage time and the trend was similar to CTC detected in the ferrofluid-selection assay, thereby validating the 100 μ l assay. The recovery of CTC from the ferrofluid-selection assay was about 90% based on the CTC in 100 μ l of blood. It was also known that MFI (Mean Fluorescence Intensity which correlates the EpCAM density) of CTC from this patient was high (MFI=300), and all EpCAM positive cells are cytokeratin positive. However, the recoveries of CTC from some other clinical samples have been as low as 20%. There may be several factors that contribute for a lower recovery, such as EpCAM positive/cytokeratin negative cells, cytokeratin dim cells, and mucin on the cell surface inhibiting the ability of ferrofluid to bind cells.

40 [0089] The assay described herein was performed on patients at two times. Response was measured by bi-dimensional imaging of the lesion. The Ratio (Ratio = Obvious CTC / Total CTC) is similar to the Response Index described earlier, and can be used as a numeric indicator of treatment success. The results are summarized in Table III. Ratios near 1.0 indicate the Total CTC are Obvious CTC, and ratios near 0.0 indicate more Suspect CTC or debris. Progressive indicates the lesion increasing in size, Partial Response indicates a response to treatment where the Ratio is relatively low, and Stabilized indicates no change, or reduction in lesion size. A positive Change indicates an increase in the number of Intact CTC, corresponding to the progression of the disease. A negative Change indicates a decrease in the number of Intact CTC, or a possible increase in the number of Suspect CTC and/or debris, corresponding to a response to treatment.

45 [0090] These results show the importance of including Suspect CTC and debris when analyzing response to treatment because the numbers of Intact or Obvious CTC alone would not provide as much information. Furthermore, such indicators are useful for short-term monitoring of treatments and therapies, or longer term monitoring for remission and/or relapse.

EP 1 425 294 B1

Table III. Obvious CTC and Suspect CTC corresponding to treatment response

Response	Ratio1	Ratio2	Change
Progressive	0.3	0.0	-0.3
	0.0	0.0	0.0
	0.5	0.6	0.1
	0.9	1.0	0.1
	0.3	0.5	0.2
	0.4	0.7	0.3
	0.0	0.4	0.4
	0.5	0.9	0.4
	0.0	0.5	0.5
	0.0	0.6	0.6
Partial Response	1.0	0.0	-1.0
	0.4	0.0	-0.4
	0.3	0.0	-0.3
	0.5	0.2	-0.3
	0.4	0.3	-0.1
	0.0	0.0	0.0
	0.3	1.0	0.7
Stabilized	1.0	0.0	-1.0
	0.5	0.0	-0.5
	1.0	0.7	-0.3
	0.3	0.0	-0.3
	0.4	0.3	-0.1
	0.1	0.0	-0.1
	0.2	0.1	-0.1
	0.6	0.5	-0.1
	0.9	0.8	-0.1

0.0	0.0	0.0
0.6	0.7	0.1
0.6	0.8	0.2
0.0	1.0	1.0

[0091] Enumeration of tumor cell debris may prove more significant in cancer diagnostics and therapeutics than detection of large proliferative cell clusters. Since debris particles in the size range, probably about 1-3 μ m (the size of platelets), have been observed to be present in much larger amounts than intact cells, they may constitute a separate, independent, and possibly more sensitive marker than intact tumor cells. The presence of damaged CTC may be particularly relevant in detecting early-stage cancer, when the immune system is intact and most active. Similarly, dramatic increases in debris during therapy may suggest breakdown of both circulating and tissue tumor cells (i.e. therapeutic effectiveness), paralleling the massive release of cellular components like calcium observed during tumor disintegration. Like soluble tumor markers, such debris may be detectable in blood without enrichment, or with minimal enrichment in the buffy coat layer and constitute an alternative, and potentially simpler diagnostic tool than intact cell enrichment/analysis. Since morphology is lost in CTC debris, detection could be done by flow cytometry as long as the debris is stained for the appropriate determinants, such as cytokeratin.

[0092] As previously discussed, damaged or fragmented CTC with or without DNA are theoretically to be expected, and therefore are not undesirable events in specimens from patients undergoing effective therapy and in untreated patients with strong immune systems. The ratio or % of intact CTC to total detectable events may prove to be a more useful parameter to the clinician in assessing a patient's immune system or response to therapy. The normal immune

EP 1 425 294 B1

defenses, especially activated neutrophils, also can damage or destroy CTC as foreign species by a process called "extracellular killing" even if the CTC are larger than the neutrophils. It does not seem surprising to find only a small percentage of the shed CTC as intact cells, unless the immune system is overwhelmed in the late stages of disease or therapy is ineffective.

[0093] Hence, there are a number of methods for *in vitro* cancer detection: conclusive detection of intact circulating cells/clusters, and inferential methods like circulating tumor debris (including total and tumor-specific RNA/DNA, and conventional soluble tumor markers. However, no method by itself may be sufficiently sensitive. Lower specificity of debris detection compared to CTC morphology may be a problem in screening that could be minimized (e.g. with triple labeling), but it may be a lesser problem in monitoring. Further statistical analysis and correlations on debris data relative to intact CTC and diagnostic stage in patients compared to normals appear worthwhile in assessing the sensitivity and specificity of debris analysis.

[0094] Examples of different types of cancer that may be detected using the compositions, methods and kits of the present invention include apudoma, choristoma, branchioma, malignant carcinoid syndrome, carcinoid heart disease, carcinoma e.g., Walker, basal cell, basosquamous, Brown-Pearce, ductal, Ehrlich tumor, *in situ*, Krebs 2, merkel cell, mucinous, non-small cell lung, oat cell, papillary, scirrhous, bronchiolar, bronchogenic, squamous cell and transitional cell reticuloendotheliosis, melanoma, chondroblastoma, chondroma, chondrosarcoma, fibroma, fibrosarcoma, giant cell tumors, histiocytoma, lipoma, liposarcoma, mesothelioma, myxoma, myxosarcoma, osteoma, osteosarcoma, Ewing's sarcoma, synovioma, adenofibroma, adenolymphoma, carcinosarcoma, chordoma, mesenchymoma, mesonephroma, myosarcoma, ameloblastoma, cementoma, odontoma, teratoma, trophoblastic tumor, adenocarcinoma, adenoma, cholangioma, cholesteatoma, cylindroma, cystadenocarcinoma, cystadenoma, granulosa cell tumor, gynandroblastoma, hepatoma, hidradenoma, islet cell tumor, leydig cell tumor, papilloma, sertoli cell tumor, theca cell tumor, leiomyoma, leiomyosarcoma, myoblastoma, myoma, myosarcoma, rhabdomyoma, rhabdomyosarcoma, ependymoma, ganglioneuroma, glioma, medulloblastoma, meningioma, neurilemmoma, neuroblastoma, neuroepithelioma, neurofibroma, neuroma, paraganglioma, paraganglioma nonchromaffin, antiokeratoma, angioma sclerosing, angiomatosis, glomangioma, hemangioendothelioma, hemangioma, hemangiopericytoma, hemangiosarcoma, lymphangioma, lymphangiomyoma, lymphangiosarcoma, pinealoma, carcinosarcoma, chondrosarcoma, cystosarcoma phyllodes, fibrosarcoma, hemangiosarcoma, leiomyosarcoma, leukosarcoma, liposarcoma, lymphangiosarcoma, myosarcoma, myxosarcoma, ovarian carcinoma, rhabdomyosarcoma, sarcoma (Kaposi's, and mast-cell), neoplasms (e.g., bone, digestive system, colorectal, liver, pancreatic, pituitary, testicular, orbital, head and neck, central nervous system, acoustic, pelvic, respiratory tract, and urogenital), neurofibromatosis, and cervical dysplasia.

[0095] However, the present invention is not limited to the detection of circulating epithelial cells and/or clusters, fragments, or debris. For example, endothelial cells have been observed in the blood of patients having a myocardial infarction. Endothelial cells, myocardial cells, and virally infected cells, like epithelial cells, have cell type specific determinants that are recognized by available monoclonal antibodies. Accordingly, the methods and the kits of the invention may be adapted to detect such circulating endothelial cells. Additionally, the invention allows for the detection of bacterial cell load in the peripheral blood of patients with infectious disease, who may also be assessed using the compositions, methods and kits of the invention. It would be reasonable to expect that these rare cells will behave similarly in circulation, and that fragments and/or debris will be present in similar conditions as those described hereinabove.

Claims

1. A method of testing for the presence of circulating tumour cells (CTCs) in a blood specimen comprising:

a. contacting the blood specimen with a formaldehyde donor and/or polyethylene glycol at the time of blood draw;
b. preparing a magnetically-labelled sample by mixing the blood specimen with colloidal magnetic particles within the range of 90 to 150 nm coupled to a first biospecific ligand, said blood specimen comprising a mixed cell population suspected of containing intact CTCs and further comprising:

- i. cell fragments derived from CTCs, or
- ii. cellular debris derived from CTCs;

said first biospecific ligand reacting specifically with said intact CTCs, and said cell fragments or said cellular debris, to the substantial exclusion of other specimen components;

c. contacting material from said magnetically-labelled sample with at least one additional biospecific ligand which specifically labels said intact CTCs, and said cell fragments or said cellular debris, to the substantial exclusion of other specimen components;

d. analyzing resulting material for the presence of said labelled CTCs, and said labelled cell fragments or said

EP 1 425 294 B1

labelled cellular debris.

2. A method for monitoring the presence of CTCs in a blood specimen comprising:

a. contacting the blood specimen with a formaldehyde donor and/or polyethylene glycol at the time of blood draw;
b. preparing a magnetically-labelled sample by mixing the blood specimen with colloidal magnetic particles within the range of 90 to 150 nm coupled to a first biospecific ligand, said blood specimen comprising a mixed cell population containing intact CTCs and further comprising:

- i. cell fragments derived from CTCs, or
- ii. cellular debris derived from CTCs;

said first biospecific ligand reacting specifically with said intact CTCs, and said cell fragments or said cellular debris, to the substantial exclusion of other specimen components;

c. contacting material from said magnetically-labelled sample with at least one additional biospecific ligand which specifically labels said intact CTCs, and said cell fragments or said cellular debris, to the substantial exclusion of other specimen components;

d. analyzing resulting material for the presence of said labelled CTCs, and said labelled cell fragments or said labelled cellular debris, the level of said labelled intact CTCs, said labelled cell fragments, and said labelled cellular debris indicating the level of CTCs in the blood specimen.

3. A method according to any preceding claim, wherein after the step of preparing said magnetically-labelled sample, said sample is subjected to a high gradient magnetic field to produce a separated magnetically-labelled fraction which is enriched for said intact CTCs, and said cell fragments or said cellular debris.

4. A method according to any preceding claim, wherein said analysis is selected from the group consisting of: multi-parameter flow cytometry, immunofluorescent microscopy, laser scanning cytometry, bright field base image analysis, capillary volumetry, spectral imaging analysis, manual cell analysis, and automated cell analysis.

5. A method according to any preceding claim, wherein said analysis further comprises classifying cell fragments or said cellular debris based on their origin as caused by apoptosis or necrosis.

6. A method according to claim 5, wherein said analysis further comprises classifying cell fragments or said cellular debris based on their origin as caused by mechanical damage, drug-induced damage or immunological damage.

7. A method according to claim 5, wherein said classification is based on at least one of the group consisting of: morphologic analysis and epitopic analysis.

8. A method according to any preceding claim, wherein said specimen contains clusters of CTCs in addition to intact CTCs.

Patentansprüche

1. Testverfahren bezüglich des Vorliegens von zirkulierenden Tumorzellen (CTC) in einer Blutprobe, das Folgendes aufweist:

a. Kontaktieren der Blutprobe mit einem Formaldehydspender und/oder Polyethylenglykol zum Zeitpunkt der Blutabnahme;

b. Vorbereiten einer magnetisch markierten Probe durch Vermischen der Blutprobe mit kolloidalen magnetischen Partikeln im Bereich von 90 nm bis 150 nm, die mit einem ersten biospezifischen Liganden verbunden sind, wobei die Blutprobe eine gemischte Zellpopulation mit dem Verdacht auf intakte CTC-Zellen und weiter Folgendes aufweist:

- i. von CTC-Zellen abgeleitete Zellfragmente; oder
- ii. von CTC-Zellen abgeleitete Zelltrümmer;

wobei der erste biospezifische Ligand spezifisch mit den intakten CTC-Zellen und den Zellfragmenten oder den Zelltrümmern reagiert, wobei andere Komponenten der Probe im Wesentlichen ausgeschlossen sind;

EP 1 425 294 B1

- c. Kontaktieren von Material aus der magnetisch markierten Probe mit mindestens einem zusätzlichen biospezifischen Liganden, der die intakten CTC-Zellen und die Zellfragmente oder die Zelltrümmer spezifisch markiert, wobei andere Komponenten der Probe im Wesentlichen ausgeschlossen sind;
- d. Analyse des sich ergebenden Materials auf das Vorliegen der markierten CTC-Zellen und der markierten Zellfragmente oder der markierten Zelltrümmer.

2. Verfahren zur Überwachung des Vorliegens von CTC-Zellen in einer Blutprobe, das Folgendes aufweist:

- a. Kontaktieren der Blutprobe mit einem Formaldehydspender und/oder Polyethylenglykol zum Zeitpunkt der Blutabnahme;
- b. Vorbereiten einer magnetisch markierten Probe durch Vermischen der Blutprobe mit kolloidalen magnetischen Partikeln im Bereich von 90 nm bis 150 nm, die mit einem ersten biospezifischen Liganden verbunden sind, wobei die Blutprobe eine gemischte, intakte CTC-Zellen enthaltende Zellpopulation und weiter Folgendes aufweist:

- i. von CTC-Zellen abgeleitete Zellfragmente; oder
- ii. von CTC-Zellen abgeleitete Zelltrümmer;
- wobei der erste biospezifische Ligand spezifisch mit den intakten CTC-Zellen und den Zellfragmenten oder den Zelltrümmern reagiert, wobei andere Komponenten der Probe im Wesentlichen ausgeschlossen sind;

- c. Kontaktieren von Material aus der magnetisch markierten Probe mit mindestens einem zusätzlichen biospezifischen Liganden, der die intakten CTC-Zellen und die Zellfragmente oder die Zelltrümmer spezifisch markiert, wobei andere Komponenten der Probe im Wesentlichen ausgeschlossen sind;
- d. Analyse des sich ergebenden Materials auf das Vorliegen der markierten CTC-Zellen und der markierten Zellfragmente oder der markierten Zelltrümmer, wobei der Pegel der markierten CTC-Zellen, der markierten Zellfragmente oder der markierten Zelltrümmer den Pegel von CTC-Zellen in der Blutprobe angibt.

3. Verfahren nach einem der vorstehenden Ansprüche, wobei nach dem Schritt der Vorbereitung der magnetisch markierten Probe die Probe einem hochgradienten Magnetfeld ausgesetzt wird, um einen separaten magnetisch markierten Bruchteil zu erzeugen, der mit den intakten CTC-Zellen und den Zellfragmenten oder den Zelltrümmern angereichert ist.

4. Verfahren nach einem der vorstehenden Ansprüche, wobei die Analyse aus der Gruppe ausgewählt wird, die Folgendes aufweist: Multiparameter-Flusszytometrie, Laserscanning-Flusszytometrie, Hellfeldbasis-Bildanalyse, Kapillarovolumetrie, Spektralbildanalyse, manuelle Zellanalyse und automatisierte Zellanalyse.

5. Verfahren nach einem der vorstehenden Ansprüche, wobei die Analyse weiter die Klassifizierung der Zellfragmente oder der Zelltrümmer auf der Grundlage ihres Ursprungs hinsichtlich der Entstehung durch Apoptose oder Nekrose aufweist.

6. Verfahren nach Anspruch 5, wobei die Analyse weiter die Klassifizierung der Zellfragmente oder der Zelltrümmer auf der Grundlage ihres Ursprungs hinsichtlich der Entstehung durch mechanische Beschädigung, drogeninduzierte Schäden oder immunologische Schäden aufweist.

7. Verfahren nach Anspruch 5, wobei die Klassifizierung auf mindestens einem Element der Gruppe basiert, die Folgendes aufweist: morphologische Analyse und epitopische Analyse.

8. Verfahren nach einem der vorstehenden Ansprüche, wobei die Probe zusätzlich zu intakten CTC-Zellen CTC-Cluster aufweist.

Revendications

1. Méthode de test pour déterminer la présence de cellules tumorales circulantes (CTC) dans un échantillon de sang comprenant :

- a. la mise en contact de l'échantillon de sang avec un donneur de formaldéhyde et/ou du polyéthylène glycol au moment de la prise de sang ;

EP 1 425 294 B1

b. la préparation d'un échantillon marqué magnétiquement en mélangeant l'échantillon de sang à des particules magnétiques colloïdales dans la plage de 90 à 150 nm couplées à un premier ligand biospécifique, ledit échantillon de sang comprenant une population cellulaire mixte suspectée de contenir des CTC intactes et comprenant en outre :

- i. des fragments cellulaires dérivés des CTC, ou
- ii. des débris cellulaires dérivés des CTC ;

ledit premier ligand biospécifique réagissant spécifiquement avec lesdites CTC intactes et lesdits fragments cellulaires ou lesdits débris cellulaires, à l'exclusion substantielle d'autres composants de l'échantillon ;

c. la mise en contact de matériau provenant dudit échantillon marqué magnétiquement avec au moins un ligand biospécifique supplémentaire qui marque spécifiquement lesdites CTC intactes, et lesdits fragments cellulaires ou lesdits débris cellulaires, à l'exclusion substantielle d'autres composants de l'échantillon ;

d. l'analyse du matériau obtenu pour déterminer la présence desdites CTC marquées et desdits fragments cellulaires marqués ou desdits débris cellulaires marqués.

2. Méthode pour contrôler la présence de CTC dans un échantillon de sang comprenant :

a. la mise en contact de l'échantillon de sang avec un donneur de formaldéhyde et/ou du polyéthylène glycol au moment de la prise de sang ;

b. la préparation d'un échantillon marqué magnétiquement en mélangeant l'échantillon de sang à des particules magnétiques colloïdales dans la plage de 90 à 150 nm couplées à un premier ligand biospécifique, ledit échantillon de sang comprenant une population cellulaire mixte contenant des CTC intactes et comprenant en outre :

- i. des fragments cellulaires dérivés des CTC, ou
- ii. des débris cellulaires dérivés des CTC ;

ledit premier ligand biospécifique réagissant spécifiquement avec lesdites CTC intactes et lesdits fragments cellulaires ou lesdits débris cellulaires, à l'exclusion substantielle d'autres composants de l'échantillon ;

c. la mise en contact de matériau provenant dudit échantillon marqué magnétiquement avec au moins un ligand biospécifique supplémentaire qui marque spécifiquement lesdites CTC intactes, et lesdits fragments cellulaires ou lesdits débris cellulaires, à l'exclusion substantielle d'autres composants de l'échantillon ;

d. l'analyse du matériau obtenu pour déterminer la présence desdites CTC marquées et desdits fragments cellulaires marqués ou desdits débris cellulaires marqués, le niveau desdites CTC intactes marquées, lesdits fragments cellulaires marqués et lesdits débris cellulaires marqués indiquant le niveau de CTC dans l'échantillon de sang.

3. Méthode selon l'une quelconque des revendications précédentes, dans laquelle, après l'étape de préparation dudit échantillon marqué magnétiquement, on soumet ledit échantillon à un champ magnétique à gradient élevé pour produire une fraction distincte marquée magnétiquement qui est enrichie en lesdites CTC intactes et en lesdits fragments cellulaires ou en lesdits débris cellulaires.

4. Méthode selon l'une quelconque des revendications précédentes, dans laquelle ladite analyse est choisie dans le groupe constitué par la cytométrie en flux multiparamètres, la microscopie par immunofluorescence, la cytométrie à balayage laser, l'analyse d'images de base en champ clair, la volumétrie capillaire, l'analyse par imagerie spectrale, l'analyse cellulaire manuelle et l'analyse cellulaire automatique.

5. Méthode selon l'une quelconque des revendications précédentes, dans laquelle ladite analyse comprend en outre la classification des fragments cellulaires ou desdits débris cellulaires en se basant sur leur origine, par exemple causée par une apoptose ou une nécrose.

6. Méthode selon la revendication 5, dans laquelle ladite analyse comprend en outre la classification des fragments cellulaires ou desdits débris cellulaires en se basant sur leur origine, par exemple causée par des dégâts mécaniques, des dégâts induits par des médicaments ou des dégâts immunologiques.

7. Méthode selon la revendication 5, dans laquelle ladite classification se fonde sur au moins l'un des groupes constitués par l'analyse morphologique et l'analyse épitopique.

EP 1 425 294 B1

8. Méthode selon l'une quelconque des revendications précédentes, dans laquelle ledit échantillon contient des grappes de CTC en plus de CTC intactes.

5

10

15

20

25

30


35


40

45

50


55

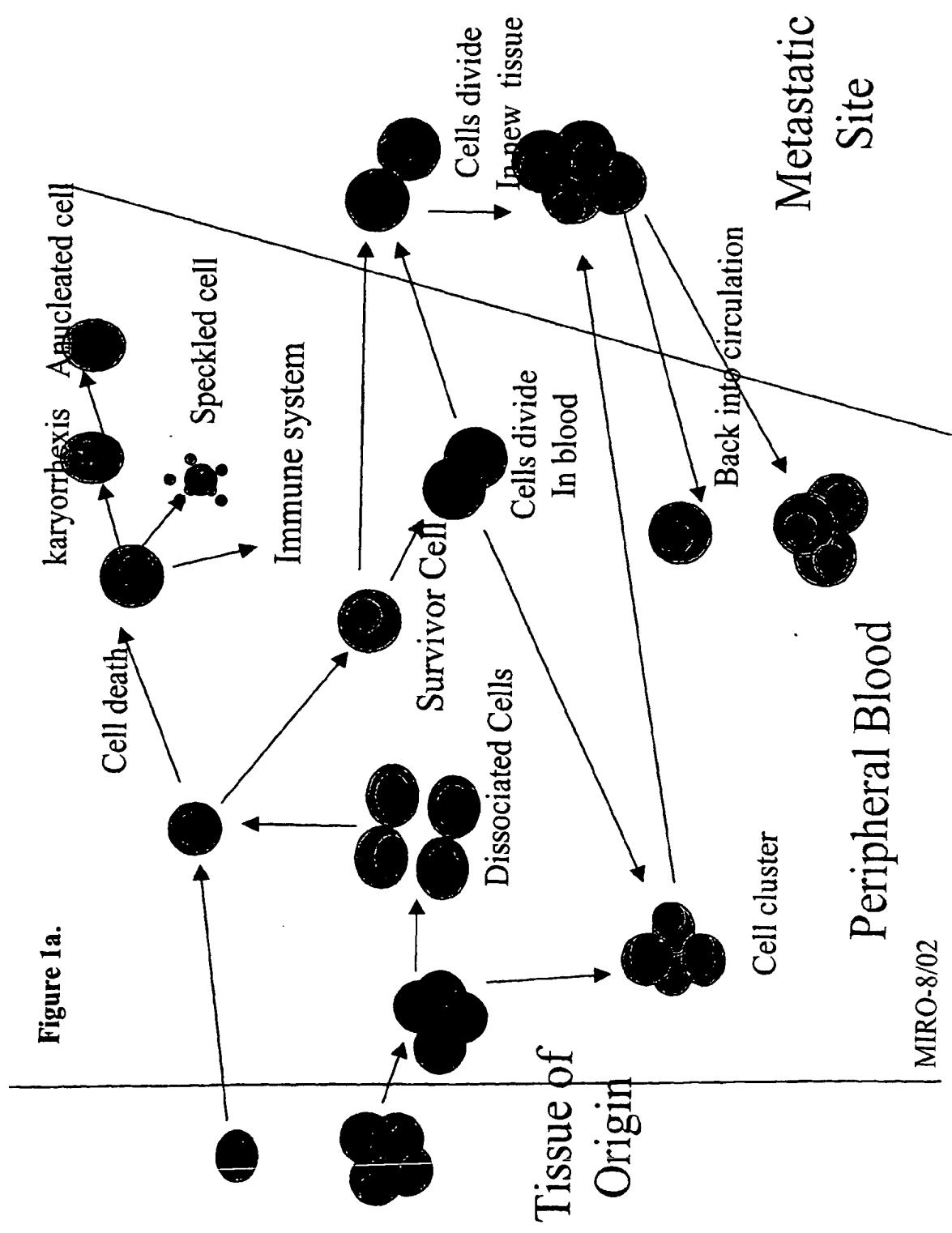

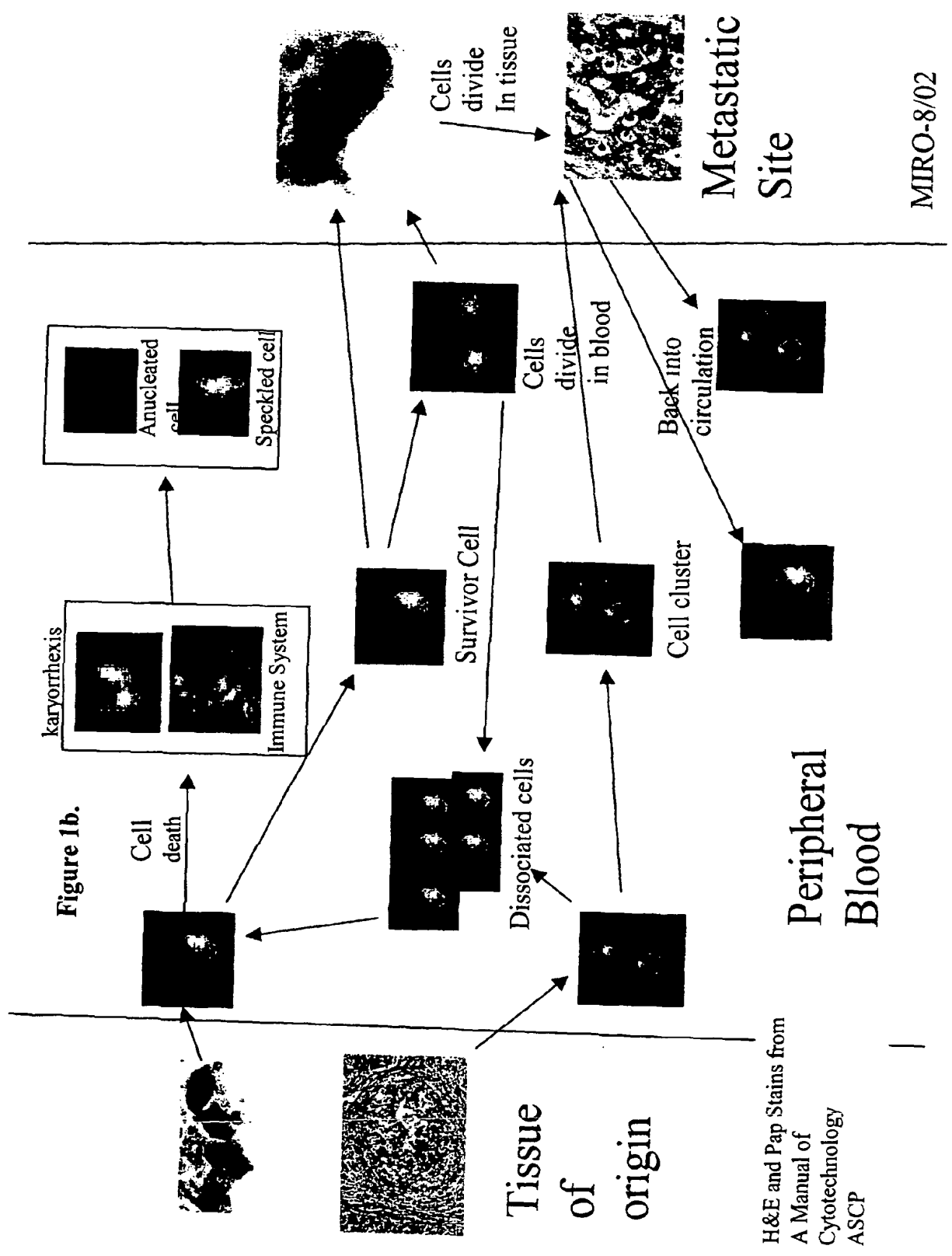


EP 1 425 294 B1

Figure 2

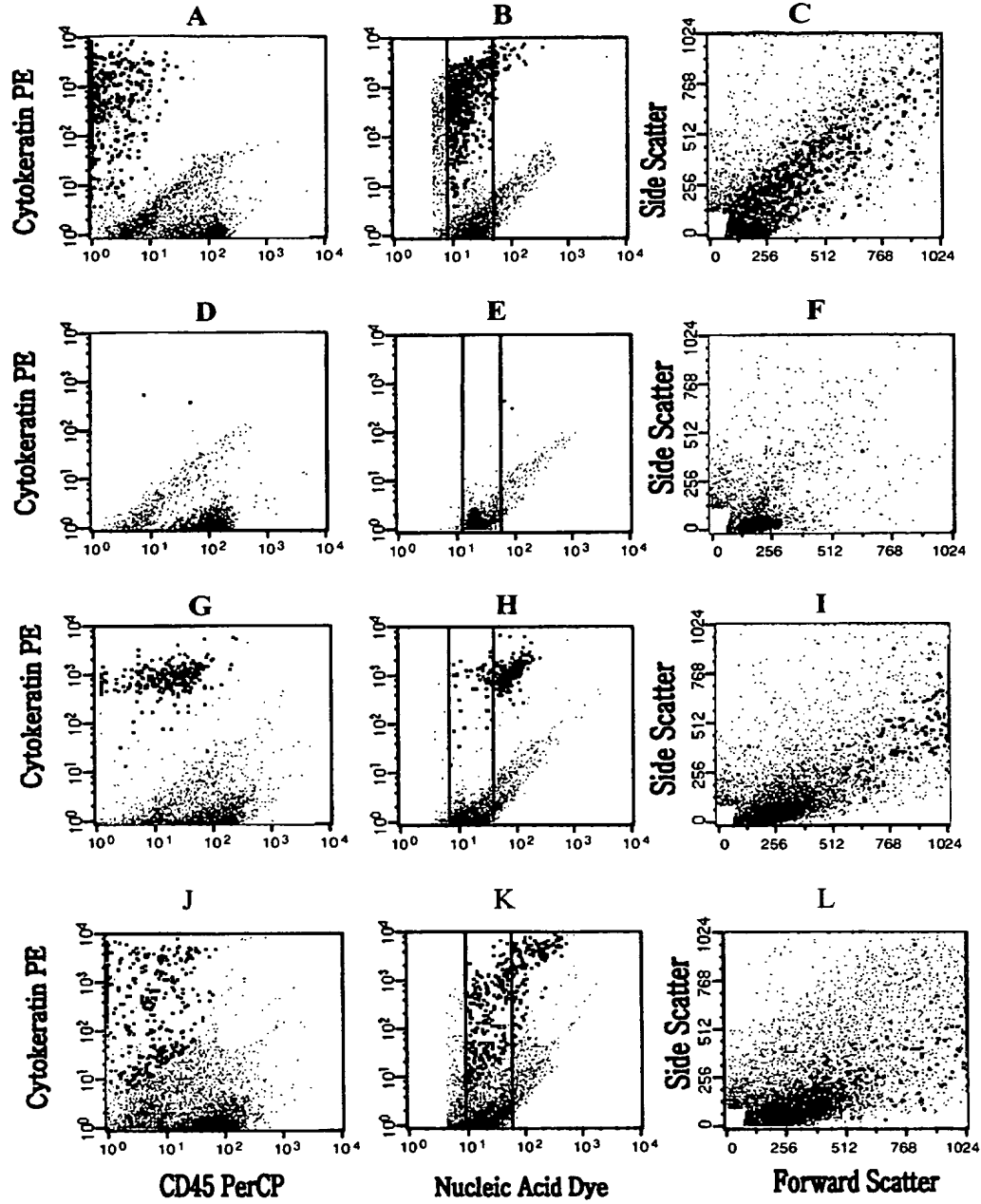

EP 1 425 294 B1

Figure 3

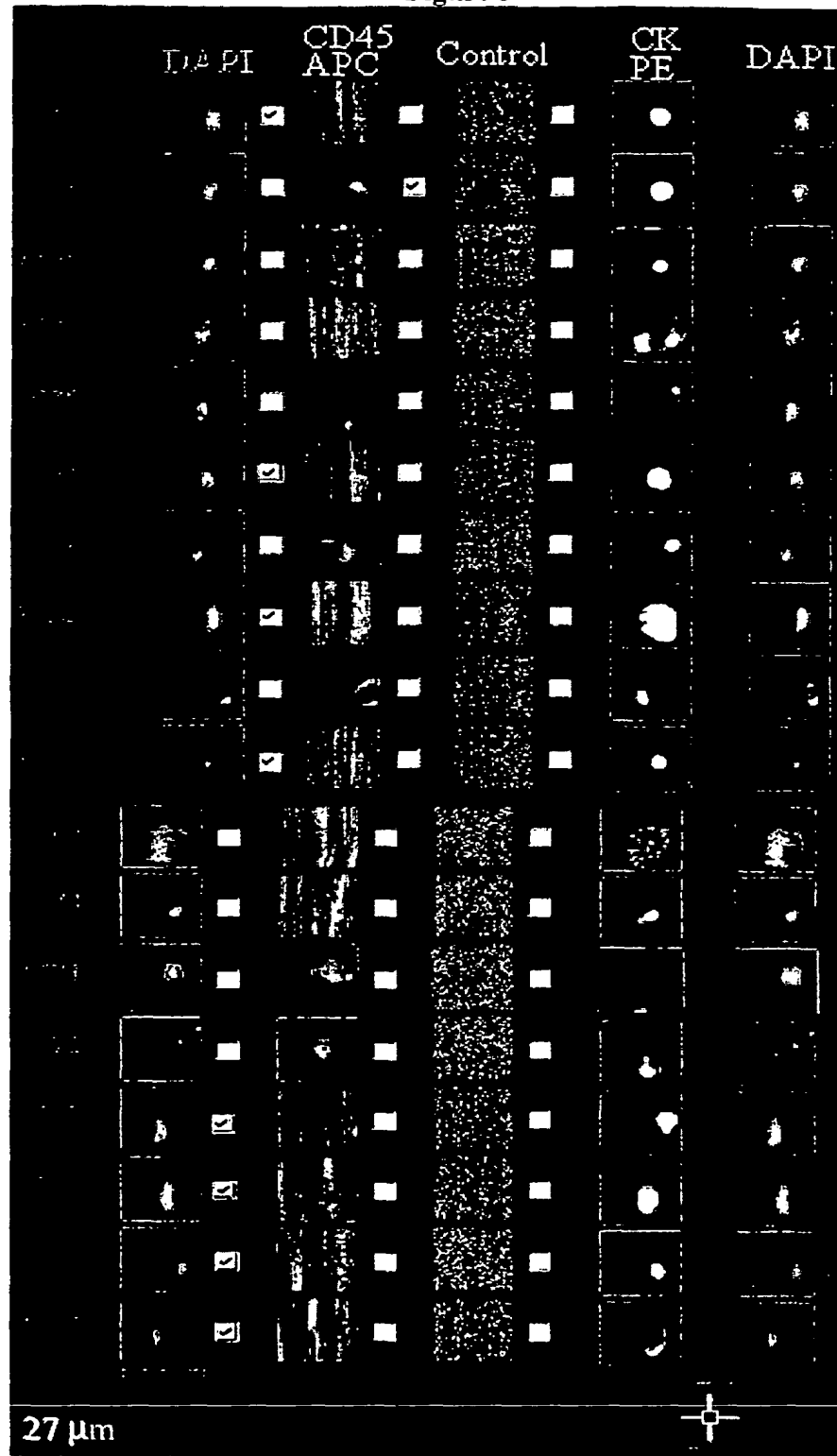

EP 1 425 294 B1

Figure 4



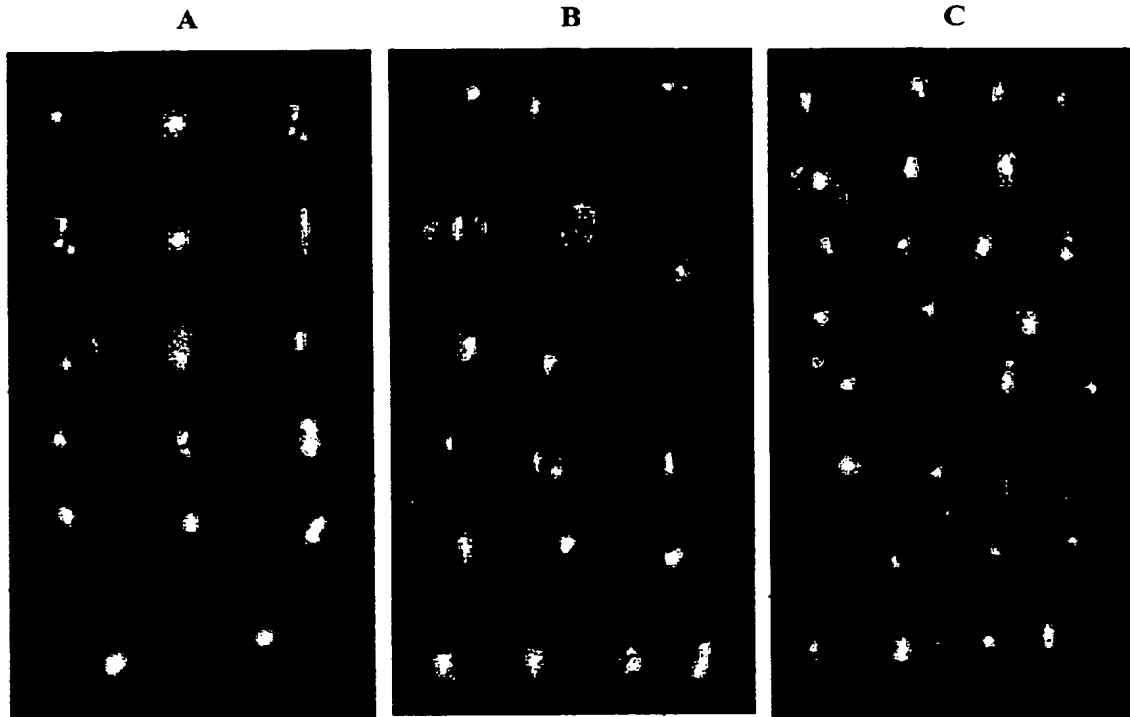


Figure 5

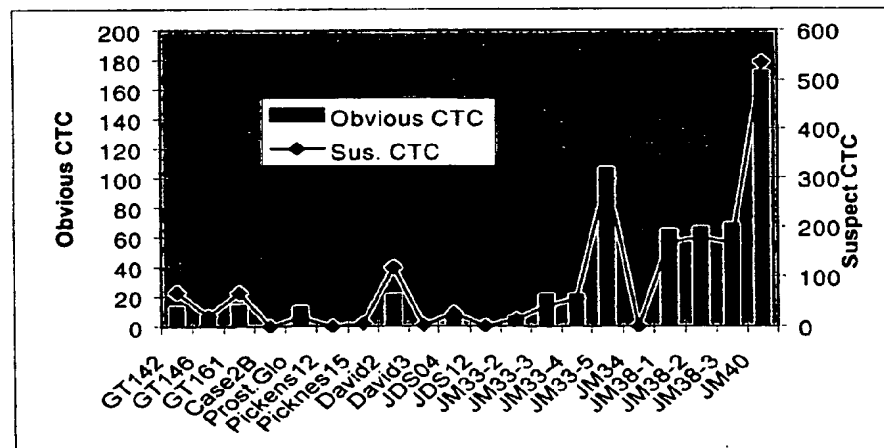

EP 1 425 294 B1

Figure 6

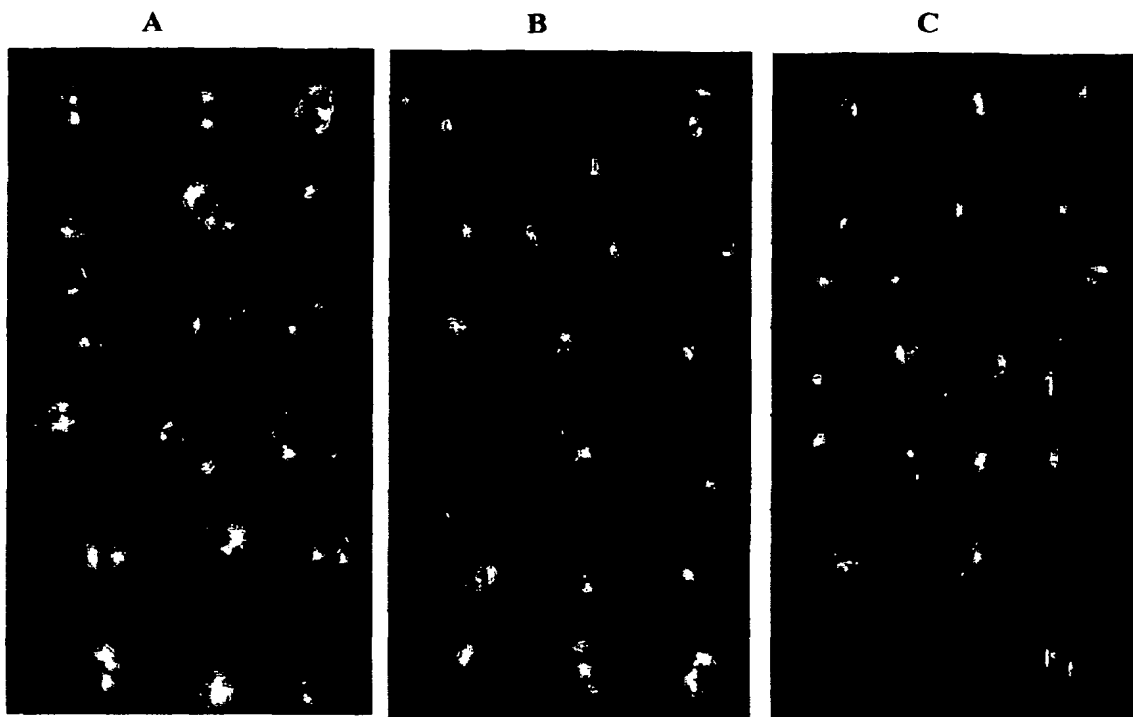

EP 1 425 294 B1

REFERENCES CITED IN THE DESCRIPTION

This list of references cited by the applicant is for the reader's convenience only. It does not form part of the European patent document. Even though great care has been taken in compiling the references, errors or omissions cannot be excluded and the EPO disclaims all liability in this regard.

Patent documents cited in the description

- US 6197523 B, Rimm [0008]
- US 6365362 B [0009] [0069]
- WO 0220825 A [0010]
- WO 0047998 A [0011] [0013] [0013]
- US 6190870 B [0013]
- WO 00479 A [0013]
- US 20010024802 A [0014]
- WO 9941613 A [0023]
- US 10081996 B [0029]
- US 3970518 A [0032]
- US 4018886 A [0032]
- US 4230685 A [0032]
- US 4267234 A [0032]
- US 4452773 A [0032]
- US 4554088 A [0032]
- US 4659678 A [0032]
- WO 03018757 A [0053]
- US 20020172987 A1 [0069]
- US 20020109838 A1 [0070]

Non-patent literature cited in the description

- **BUTLER, TP ; GULLINO PM.** *Cancer Research*, 1975, vol. 35, 512-516 [0003]
- **B BRANDT et al.** Isolation of prostate-derived single cells and cell clusters from human peripheral blood. *Cancer Research*, 1996, vol. 56, 4556-4561 [0018]
- Manual of Cytology. American Society of Clinical Pathologists Press, 1983 [0059]